

When you need a choice.

In a world where no two patients, environments or
Reusable

needs are alike, you need a choice in reusable sensors.

Our line of Nellcor® reusable sensors offers you a

comprehensive choice to effectively meet your

patient monitoring needs. From the Dura-Y®

sensor – our premium-performance reusable –

to the convenient clip design of the Durasensor,®

you're assured of the quality and durability

Mallinckrodt builds into all of its sensors.




MALLINCKRODT

NELLCOR REUSABLE SENSORS

"If the cable is gray, don't throw it away!"

In order to avoid the possibility of mistakenly throwing away a reusable sensor, Mallinckrodt manufactures all its reusable pulse oximetry sensors with gray cables?"

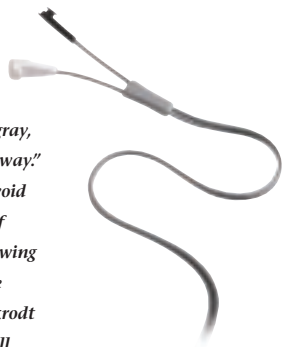

When you need a reusable sensor

You'll find that our reusable sensors are a perfect match to your monitoring needs when patients have normal perfusion and pulse, and the risk of cross-contamination is low. Reusable sensors are ideal for spot checks and short-term monitoring, although you must check the site every four hours during continuous monitoring, and relocate the sensor as necessary.

And to help standardize your sensor selection, Nellcor reusable sensors are compatible with many leading manufacturers of patient monitoring systems**

Nellcor reusable sensors – the right choice for all your reusable sensor needs.

Product Features

Dura-Y Multisite Sensor

Model D-YS

- Shipped nonsterile, 1 D-YS and a wrap assortment pack containing 40 wraps
- Built-in shielding protects the signal from electronic noise
- Bright, high quality LEDs maximize tracking capabilities
- Easy-to-use button design secures adhesive wrap to sensor
- Black-and-white sensor head design optimizes signal and shields detector from ambient light
- For patients weighing more than 1 kg
- Preferred application sites:
 - Adult (>40 kg) – index finger
 - Pediatric (15-40 kg) – index finger
 - Infant (3-15 kg) – great toe
 - Neonate (1-3 kg) – ball of the foot



Durasensor Finger Clip Sensor

Model DS-100A

- Shipped nonsterile, 1 DS-100A
- For patients weighing more than 40 kg
- Easy-to-use clip design
- Preferred application site is index finger



Oxiband® Sensor

Model OXI-A/N

- Shipped nonsterile, 1 OXI-A/N and 1 box of 50 adhesive wraps
 - Easy-to-use "wraparound" sensor design
- ADULT APPLICATION (>40 kg)
- Preferred application site is index finger
- NEONATAL APPLICATION (<3 kg)
- Preferred application site is ball of the foot



Dura-Y Ear Clip

Model D-YSE

- Shipped nonsterile, 1 D-YSE
- Accessory for Dura-Y sensor, sold separately
- Easy-to-use button design secures clip to sensor
- For patients weighing more than 30 kg
- Preferred application sites are ear lobe and pinna



Oxiband Sensor

Model OXI-P/I

- Shipped nonsterile, 1 OXI-P/I and 1 box of 50 adhesive wraps
 - Easy-to-use "wraparound" sensor design
- PEDIATRIC APPLICATION (15-40 kg)
- Preferred application site is index finger
- INFANT APPLICATION (3-15 kg)
- Preferred application site is great toe



Pedi Check™ Pediatric Spot-Check Clip

Model D-YSPD

- Shipped nonsterile, 1 D-YSPD
- Accessory for Dura-Y sensor, sold separately
- For patients weighing 3 to 40 kg
- Easy-to-use clip design for attended spot-check monitoring not to exceed 20 minutes
- Preferred application site is index finger



RS-10 Reflectance Sensor

Model RS-10

- Shipped nonsterile, case of 6 with 6 headbands
- For patients weighing more than 40 kg
- Limited reuse
- Preferred application site is on the forehead above the brow or below the hairline



For patients with compromised peripheral circulation or long-term monitoring needs, please consult your local Mallinckrodt representative about Oxisensor® II sensors.

*The RS-10 Reflectance Sensor is a "limited reuse" sensor and is manufactured with a white cable.

** Check individual manufacturers for compatibility of particular instruments and sensor models.

© 1999 Mallinckrodt Inc. All rights reserved. 02771-1099

MALLINCKRODT

Mallinckrodt Inc.
675 McDonnell Boulevard
P.O. Box 5840
St. Louis, MO 63134 USA
Tel 314.654.2000
Toll Free 1.800.635.5267
www.mallinckrodt.com

Mallinckrodt
Europe BV
Hambakenwetering 1
5231 DD 's-Hertogenbosch
The Netherlands
Tel +31.73.6485200

Mallinckrodt
Asia Pacific Pte Ltd.
83 Clemenceau Avenue
#16-03 UE Square
Singapore 239920
Tel +65.738.4333