


Healthcare is changing.
Are you prepared to change with it?

IQ
stress



A woman is shown from the waist up, running on a treadmill. She is wearing a dark tank top and has several white ECG electrodes attached to her chest and shoulders. A blue monitoring device is clipped to the treadmill's handrail, with several blue cables connected to it. The background is dark, and the overall lighting is dramatic, highlighting the woman's form and the equipment.

Coronary Heart Disease (CHD) is the single largest cause of death in the United States. In fact, almost one in three adults suffer from some form of cardiovascular disease. Proper screening and care for cardiac diseases is an essential part of any patient care plan. By adding ECG and stress testing to your diagnostic assessment, your patients can be assured of your commitment to providing the care they need, whenever they need it.

DISCOUNT CARDIOLOGY
discountcardiology.com

Simplified stress testing for complete patient care

IQ
stress

The IQstress™ system from Midmark is an ideal component of any complete patient care strategy. It combines sophisticated cardiac testing with easy-to-use features aimed at increasing the time you spend with your patients.





The lightweight ECG module can be worn comfortably by the patient or easily mounted on the treadmill.

“Thanks to the IQstress System for increasing efficiencies within our practice. Since implementation of the product, we have increased our procedure throughput by 17%, all while producing high quality diagnostic stress studies for our patients.”



Michelle Roesler, Director of Ancillary Services,
Buffalo Medical Group, PC, Yorba Linda, CA

The stress testing workflow goes digital


IQstress™ Digital Stress Testing System

The IQstress system from Midmark is an affordable, easy-to-use system that combines sophisticated cardiac testing with straightforward features aimed at increasing the time you spend with your patients. The IQstress system includes a lightweight, digital ECG, integrated stress analysis software, stress cart and a medical grade treadmill. Features include:


- Everything is just a click away with a user-friendly, one-level menu
- Customizable user profiles make procedure setup fast and easy
- Quick access to customized test protocols, measurement parameters, display formats and reports
- Seven protocols pre-programmed to begin testing immediately, including the Bruce and Modified Bruce
- Easily review up to four stress tests for side-by-side comparison of a patient's progress
- Summarize and review test results with detailed rhythm strips and averaged beats immediately after the test is completed, either locally or from any network
- The medical grade treadmill is reliable, quiet and runs smoothly for assured patient comfort

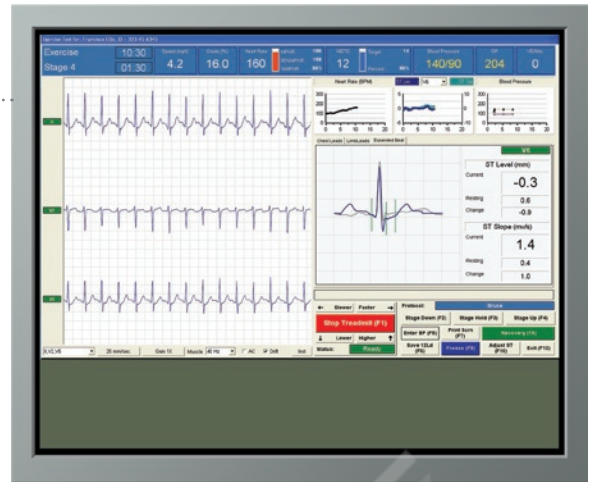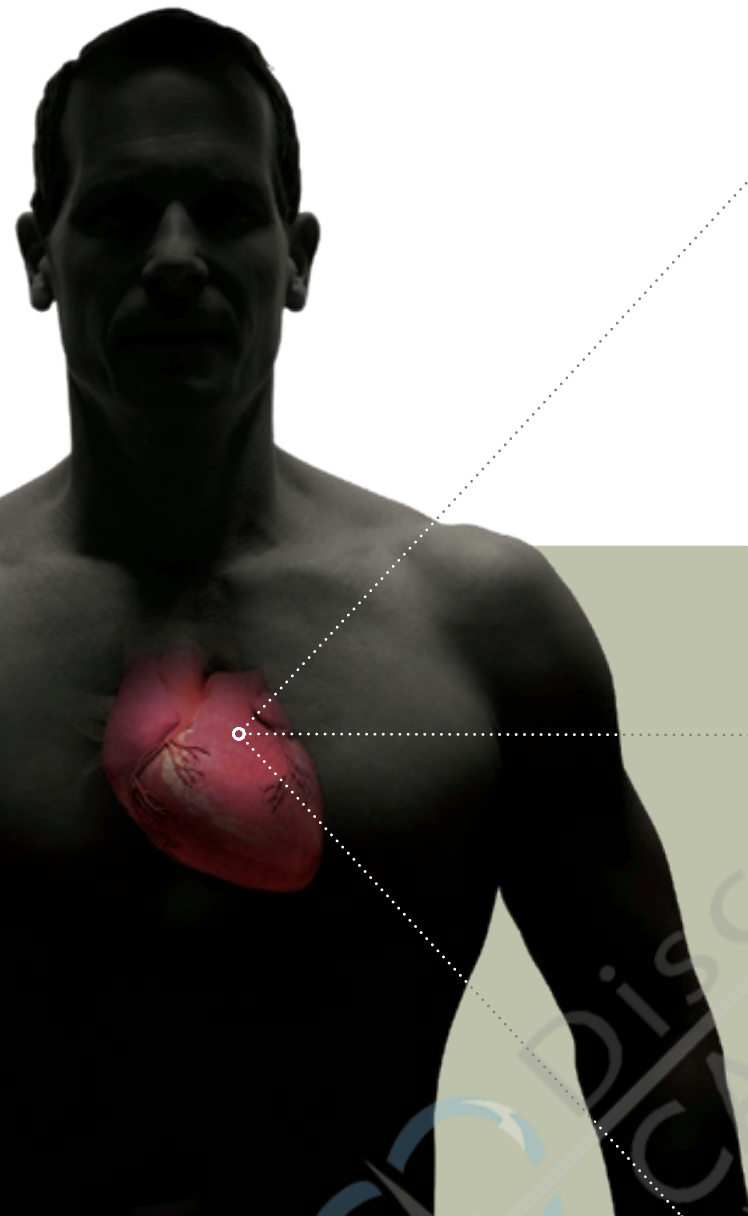

Live Monitoring in Exercise Phase

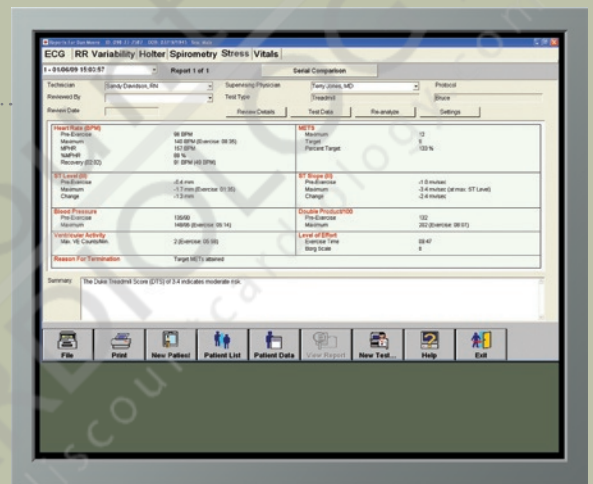

Stress Test Result Summary

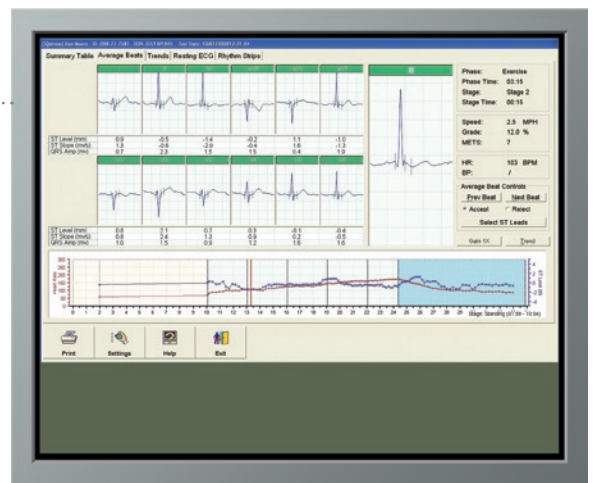

Average Beats Trend and Review

IQmanager™ Software

Each Midmark IQstress system includes the IQmanager software, a simple, cost effective way to capture and store patient test results right on your existing computer system. It eliminates the need for thermal paper and reduces the copying, scanning and shredding of paper records.

- Built in patient management software to review, edit and store patient information from any computer on your network
- Acquires, organizes and presents patient data for comprehensive view
- Simple to use with one click email and fax options

- Save time by eliminating need to re-enter patient data with easy-to-use database management software
- Email or fax test reports directly from your computer or archive them in CAR, PDF and JPG formats.

EHR/EMR Integration

Integrate test results and patient data directly into your EHR/EMR solution.

- Seamless integration allows you to perform tests within the EHR or EMR and store the information in the patient's chart
- Works with more than 80 of today's top EHR/EMR solutions

Make your practice more successful today with an IQ Digital Diagnostic System

Productivity gains and improvements in patient care are just a few of the benefits that can result from implementing the IQ Digital Diagnostic System. You'll have the ability to provide more streamlined diagnostic procedures for your patients and a more efficient workflow for your practice. Call us today and we'll put together a customized suite of solutions to help your practice move toward the future of healthcare.



IQecg™



IQspiro™



IQholter™, EX, EP



IQmanager™



IQstress™

Ordering Information

- 4-000-0332 IQstress™ w/Treadmill & Cart
- 4-000-0334 IQstress™ w/Treadmill, Cart, Computer, Printer, UPS
- 2-100-0095 IQmark Stress Prep Kits (box of 5)
- 4-100-1200 IQmanager™ Software - Patient data management of ECG, Spirometry, Holter and Stress
- 3-100-0202 IQstress™ Patient Cable - blue
- 3-370-0010 IQstress™ Pouch & Belts

For information, pricing or to place an order, contact us at:
1-800-624-8950 or
1-310-516-6050 or
midmark.com

Computer analysis and interpretation is not a substitute for physician interpretation and, therefore, should be examined with respect to the patient's overall clinical condition.

Midmark Diagnostics Group is an ISO 13485 Certified Company.
Midmark Diagnostics Group
1125 W. 190th Street, Gardena, CA 90248
midmark.com
1-800-624-8950

© 2009 Midmark Corporation
60 Vista Drive P.O. Box 286 Versailles, Ohio 45380-0286
Products subject to improvement changes without notice
Litho in U.S.A. 007-0744-00 Rev. A1 (2/09)



Because we care.