

P10

Digital
Ultrasonic
Diagnostic
Imaging
System



P10



Characterized by a compact lightweight design with user-centric workflow, the P10 delivers on the need for highly portable easy-to-use ultrasound.

Features

The P10 is a portable Diagnostic Ultrasound System, which applies advanced technologies such as Phased Inversion Harmonic Compound Imaging (eHCI), Multi-Beam-Forming (mBeam), Speckle Resistance Imaging (eSRI), and Spatial Compounding Imaging, etc. Various image parameter adjustments, 12.1 inch LCD and diverse probes are configured to provide clear and stable images.

- Multi-language displays
- 256 gray scale
- Two probe sockets

Display modes:

B, B+B, 4B, B+M, M, B+PW, B+Color, B+Color+PW, B+PDI/DPDI, B+PDI/DPDI +PW.

Generic Measurements and Calculations:

- B/Color Mode: Distance, Cir/Area (Ellipse/Trace), Volume (2-Axis/3-Axis), Ratio, % Stenosis, Angle, and Histogram.
- M Mode: Distance, Time, Slope and Hear Rate (two cycles).
- PW Mode: Velocity, Heart Rate, Time, Acceleration, Resistance Index (RI), Pulsatility Index (PI) and Auto (auto trace).

Technical Specifications

Display Modes

B Mode:	Single, Dual, Quad
C Mode:	B/C (Single, Dual);
Simultaneous dual mode:	B+B/C; Triplex mode B/C/PW
PDI/DPDI Mode:	B/PDI (DPDI) (Single, Dual); Simultaneous dual mode: B+B/PDI(DPDI); Triplex mode B/PDI(DPDI)/PW
PW Mode:	B/PW (duplex, simultaneous); B+C/PW, B+PDI (DPDI) / PW (duplex); B/C/PW, B/PDI (DPDI) / PW (triplex);
M Mode:	B/M (Display layout: Up/down, Left/right, 1:1)
Image Gray Scale:	256 levels
Image Magnification:	In area Real time: ×1.2, ×1.4, ×1.6, ×2.0, ×2.4, ×3.0, ×4.0 Frozen: ×1.14, ×1.33, ×1.6, ×2.0, ×2.67, ×3.2, ×4.0
Storage:	504 MB
Cine Review:	409 frames (Color) / 1227frames (Black & White)
Depth Adjustment:	Adjustable in real time in all modes
Image Conversion:	Up/Down flip, Left/Right flip, 90° rotate, B/W invert
Language Conversion:	Chinese, English, etc. (The language options varies with language software installed.)
Focus Number:	Max. 4
Software Packages:	Abdomen, obstetric, small parts, gynecology, cardiology, urology, vascular and Pediatrics.
B mode Measurement:	Distance, circumference, area, volume, ratio, % stenosis, angle, and histogram.
M mode Measurement:	Distance, time, slope, and heart rate
D mode Measurement:	Time, heart rate, velocity, acceleration, RI, PI and Auto (auto trace)
Annotations:	Patient name, age, sex, time, date, hospital name, doctor name, comment (full-screen character editing)
Body Mark:	130 types
USB Port:	USB 1.1

Transducers

Convex array:
C352UB



Micro-convex array:
C612UB



Micro-convex array:
C6152UB



Linear array:
L552UB



Micro-convex array:
C422UB



Linear array:
L742UB



Linear array:
L1042UB




Endovaginal:
E612UB




Screen Samples



Color Doppler Sensitivity
Renal Flow

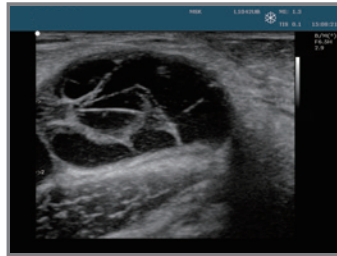


Harmonic Imaging
Baker's Cyst

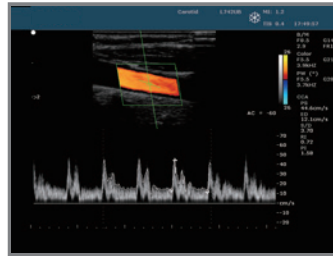

Auto-Doppler Calculations
Carotid Artery



Endovaginal Imaging
First Trimester OB

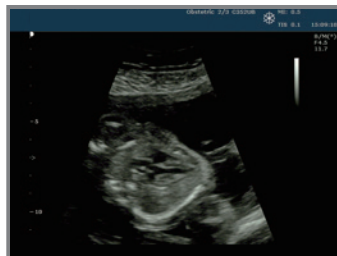

2D Spatial Resolution
Fetal Heart



Color Doppler Imaging
Second Trimester OB



Color Spatial Resolution
Hepatic Veins



MEDISONO,
3511 Silverside RD.
Wilmington, DE 19810
United States
info@medisono.com
www.medisono.com