



DIAGNOSTIC PRODUCTS

Simple Solutions, Significant Results™



DIAGNOSTICS

Blood Pressure Monitors

Medline Diagnostics

Medline is proud to present our latest line of diagnostics. We have built a reputation on providing diagnostics that are highly reliable, extremely accurate, and built to last.



Index

Blood Pressure Monitors.....	3 – 5
Sphygmomanometers.....	6 – 7
Accessories.....	8 – 9
Stethoscopes.....	10 – 13
Compass Diabetes Program.....	14 – 15
Point of Care Testing.....	16 – 17
Lancets.....	17 – 18
Sharps Containers.....	20 – 23
Patient Monitoring.....	24 – 25
Thermometers.....	26 – 30
Defibrillators.....	31
Bladder Scanner.....	32

If you have questions about any of the products in this brochure, contact your Medline representative or call us at 1-800-MEDLINE.

DIGITAL BLOOD PRESSURE MONITORS



MDS2007

Desktop Blood Pressure Monitor Features Average Mode Technology!

These premium monitors detect the appearance of an irregular heartbeat during measurement and gives a warning signal with your reading.

- Measurement includes patented Irregular Heartbeat Detector.
- Monitor also automatically stores 99 readings with date and time.
- Average mode feature takes three consecutive measurements and averages them for more accurate results.

Item #	Description	Pkg.
MDS2007	Premium Monitor. Includes 4-AA batteries, medium arm cuff and deluxe storage case.	1 ea

Pressure Monitor Accessories Replacement cuffs and adapters.

Item #	Description	Pkg.
MDS9872	Adult Cuff, (fits arms 22-30 cm, 8.58"-11.70")	1 ea
MDS9873	Large Adult Cuff, (fits arms 30-38 cm, 11.70"-14.82")	1 ea
M-IADE01	AC Adapter	1 ea



Real-time graphic display of blood pressure being measured, on-screen instructions.



MDS2006

Average Mode Technology Wrist Blood Pressure Monitor

This advanced blood pressure monitor uses a unique, patented Average Mode Technology to take three consecutive measurements then display a single, averaged result — a procedure practiced by many physicians for a more accurate result.

- Real-time graphic display of blood pressure being measured, on-screen instructions.
- 5 languages.
- 30 memories for 2 separate users with date and time.
- Average mode feature takes three consecutive measurements and averages them for more accurate results.

Item #	Description	Pkg.
MDS2006	Wrist Monitor	1 ea



DIGITAL BLOOD PRESSURE MONITORS

Blood Pressure Monitors



MDS2003

Digital readout of systolic, diastolic pressure and pulse

Wrist Blood Pressure Monitor

- One-button operation
- Automatically inflates and deflates
- Digital readout of systolic and diastolic pressure
- Digital readout of pulse
- Memory with date and time
- Operates on two AAA batteries

Item #	Description	Pkg.
MDS2003	Digital BP Wrist Monitor	1 ea

Automatic Digital Blood Pressure Unit

Our revolutionary automatic digital blood pressure units have automatic inflation and error detection thanks to the proprietary "fuzzy logic" programming. The easy one-button design provides accurate and reliable readings.

- "Fuzzy logic" circuitry allows this unit to customize inflation levels for individual patients
- Eight-bit processor provides unsurpassed accuracy
- Auto pulse readings
- Memory
- Large digital readout
- Operates on four-AA batteries



MDS2001

Item #	Description	Pkg.
MDS2001	Automatic Inflation w/Std Adult Cuff (fits arms 22-30 cm, 8-12 inches)	1 ea
MDS2001LA	Automatic Inflation w/Lrg Adult Cuff (fits arms 30 -38 cm, 12-15 inches)	1 ea
MDS9873	Large Adult Cuff (fits arms 30-38 cm, 12-15 inches)	1 ea
MDS9872	Adult Cuff	1 ea
M-IADE01	AC Adapter	1 ea

Digital Blood Pressure Unit – Manual Inflate

Manually inflated digital blood pressure units have the convenience of a large, digital readout. Also provides accurate and reliable readings.

- Manual inflation
- Auto pulse readings
- Memory
- Large digital readout
- Operates on four-AA batteries



MDS2002

Item #	Description	Pkg.
MDS2002	Manual Inflation w/Std Adult Cuff (fits arms 22-30 cm, 8-12 inches)	1 ea
MDS2002LA	Manual Inflation w/Lrg Adult Cuff (fits arms 30-38 cm, 12-15 inches)	1 ea
MDS9873	Large Adult Cuff (fits arms 30-38 cm, 12-15 inches)	1 ea
MDS9872	Adult Cuff	1 ea
M-IADE01	AC Adapter	1 ea



DIGITAL BLOOD PRESSURE MONITORS

Large, raised numbers and a luminescent (glow in the dark) coating make the Nite Shift Wall-Mount Aneroid easy to read in almost any type of lighting.

- Eight-foot coiled tubing and pivoting bracket extend service area for more convenient use.
- Inflation system and bulb assembly store easily behind gauge.
- Features Medline's nylon range finder cuff.
- Available in choice of Latex or Latex-Free

Item #	Description	Pkg.
MDS9400	Adult, Latex	1 ea
MDS9400LF	Adult, Latex-Free, Neoprene	1 ea
MDS9472	Wall mount basket (not shown)	1 ea
MDS9418	Mobile Stand	1 ea



MDS9400LF



MDS9407LF

Blood Pressure Monitors

The Mobile Aneroid maintains all of the features of its Wall-Mount counterpart, the 9400, while adding the convenience of portability and height adjustment.

- Adjusts from 29" to 48" (74 cm to 122 cm) to accommodate a wide range of viewing levels.
- Includes inflation system and bulb assembly storage basket.

Item #	Description	Pkg.
MDS9407	Adult, Latex	1 ea
MDS9407LF	Adult, Latex-Free Neoprene	1 ea



Blood Pressure Cuff Barriers

Single-use wrap for added protection from cross-contamination.

- Protection against bacteria pass through
- Will not alter blood pressure readings
- Latex-Free
- Available in 3 sizes to fit any size blood pressure cuff

Item #	Description	Pkg.
MDS916113	Small	300/cs
MDS916114	Adult	300/cs
MDS916115	Large	150/cs



Protection from Cross Contamination

Blood Pressure Cuff-Guard

Help reduce the potential of cross-contamination by using our protective, disposable blood pressure cuff covers.

- Latex-free
- Made of Tyvek® material
- Quick and easy to replace or remove the cover!

Item #	Description	Pkg.
MDS9157	Child Cuff-Guard	20/bx
MDS9158	Adult Cuff-Guard	20/bx
MDS9159	Large Adult Cuff-Guard	20/bx
Soft Fabric for Single Time Use		
MDS9167	Child Soft Cuff-Guard	20/bx
MDS9168	Adult Soft Cuff-Guard	20/bx
MDS9169	Lg. Adult Cuff Guard	20/bx

SPHYGMOMANOMETERS



MDS9114

The Compli-Mate Aneroid Sphygmomanometers are designed for many years of reliable service.

- Designer-colored cuffs offer a choice of vibrant colors.
- Comes complete with oversized matching nylon carrying case which allows the user to carry additional items.

Item #	Description	Pkg.
MDS9104	Red	1 ea
MDS9106	Navy Blue	1 ea
MDS9107	Lavender	1 ea
MDS9108	Hunter Green	1 ea
MDS9109	Orchid	1 ea
MDS9112	Teal	1 ea
MDS9113	Purple	1 ea
MDS9114	Royal Blue	1 ea

Combining the Compli-Mates Aneroid Sphygmomanometers with the very popular Sprague Rappaport Stethoscope makes this Combination Kit popular with nurses and EMT's.

- Kit comes with an aneroid sphygmomanometer and a 22" (56 cm) Sprague Rappaport stethoscope.
- Packed together in a matching nylon carrying case.

Item #	Description	Pkg.
MDS9121	Magenta	1 ea
MDS9123	Purple	1 ea
MDS9124	Royal Blue	1 ea
MDS9125	Black	1 ea
MDS9127	Lavender	1 ea
MDS9128	Hunter Green	1 ea

The Compli-Mates Dual-Head Combination Kit combines an aneroid sphygmomanometer and a dual-head stethoscope.

- Ideal for nurses and students.
- Kit comes complete with an aneroid sphygmomanometer and matching dual-head stethoscope.
- Packed together in a matching nylon carrying case.

Item #	Description	Pkg.
MDS9116	Magenta	1 ea
MDS9117	Teal	1 ea
MDS9118	Purple	1 ea
MDS9119	Royal Blue	1 ea
MDS9133	Hunter Green	1 ea
MDS9135	Burgundy	1 ea
MDS9136	Red	1 ea



MDS9117

SPHYGMOMANOMETERS

The Latex-Free Handheld Aneroid's lightweight design and portability make it a popular handheld unit.

- No-stop pin action allows a visual calibration check, ensuring accuracy.
- Nylon rangefinder cuff provides quick, convenient patient fitting.
- Latex-Free construction.
- Field adjustable for long-term accuracy.

PVC			Neoprene		
Item #	Description	Pkg.	Item #	Description	Pkg.
MDS9380	Adult	1 ea	MDS9380LF	Adult	1 ea
MDS9386	Infant	1 ea	MDS9386LF	Infant	1 ea
MDS9387	Child	1 ea	MDS9387LF	Child	1 ea
MDS9388	Large Adult	1 ea	MDS9388LF	Large Adult	1 ea



MDS9380LF

Our Best Selling Sphygmomanometer.

Our finest handheld aneroid, the Nite-Shift Premier, has a luminescent dial that facilitates reading over a wide range of lighting conditions.

- Durable zinc-alloy shell protects precision aneroid mechanism.
- Shatter resistant acrylic faceplate surrounds the luminescent dial.
- Field adjustable for long-term accuracy.

Item #	Description	Pkg.
MDS9410	Adult	1 ea
MDS9411	Child	1 ea
MDS9413	Large Adult	1 ea



MDS9410

Sphygmomanometers

FOR HOME USE

Medline's Home Blood Pressure Kits feature D-ring cuffs to simplify home use

- Both units feature an easy-to-read aneroid gauge, complete instructions, recording chart, and carrying case.
- MDS9301 features an attached stethoscope while the two-person MDS9300 features a separate single-head stethoscope.
- Ideal for home use.

Item #	Description	Pkg.
MDS9300	With detached stethoscope	1 ea
MDS9301	With attached stethoscope	1 ea

Simple and convenient



MDS9301


ACCESSORIES



MDS91420LF



MDS91410

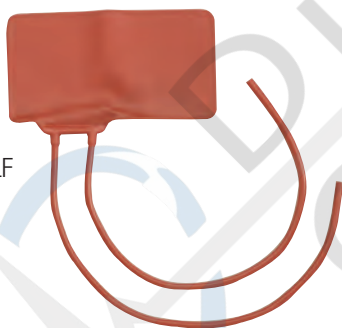

MDS9505LF



MDS9205



MDS9420

Two-Tube Inflation Bag & Nylon Range finder Cuff

Item #	Description	Size	Pkg.
Latex-Free (PVC)			
MDS91420	Adult	5 $\frac{3}{4}$ " x 22", (15 cm x 56 cm)	1 ea
MDS91421	Child	4 $\frac{1}{4}$ " x 15", (11 cm x 38 cm)	1 ea
MDS91422	Infant	2 $\frac{7}{8}$ " x 10 $\frac{1}{4}$ ", (7 cm x 26 cm)	1 ea
MDS91423	Large Adult	6 $\frac{7}{8}$ " x 26", (17 cm x 66 cm)	1 ea
MDS91424	Thigh	8 $\frac{1}{2}$ " x 30 $\frac{1}{2}$ ", (22 cm x 77 cm)	1 ea
Neoprene Latex-Free			
MDS91420LF	Adult	5 $\frac{3}{4}$ " x 22", (15 cm x 56 cm)	1 ea
MDS91421LF	Child	4 $\frac{1}{4}$ " x 15", (11 cm x 38 cm)	1 ea
MDS91422LF	Infant	2 $\frac{7}{8}$ " x 10 $\frac{1}{4}$ ", (7 cm x 26 cm)	1 ea
MDS91423LF	Large Adult	6 $\frac{7}{8}$ " x 26", (17 cm x 66 cm)	1 ea

One-Tube Inflation Bag & Nylon Range finder Cuff

Item #	Description	Size	Pkg.
Latex-Free (PVC)			
MDS91410	Adult	5 $\frac{3}{4}$ " x 22", (15 cm x 56 cm)	1 ea
MDS91411	Child	4 $\frac{1}{4}$ " x 15", (11 cm x 38 cm)	1 ea
MDS91412	Infant	2 $\frac{7}{8}$ " x 10 $\frac{1}{4}$ ", (7 cm x 26 cm)	1 ea
MDS91413	Large Adult	6 $\frac{7}{8}$ " x 26", (17 cm x 66 cm)	1 ea
MDS91414	Thigh	8 $\frac{1}{2}$ " x 30 $\frac{1}{2}$ ", (22 cm x 77 cm)	1 ea
Neoprene Latex-Free			
MDS91410LF	Adult	5 $\frac{3}{4}$ " x 22", (15 cm x 56 cm)	1 ea
MDS91411LF	Child	4 $\frac{1}{4}$ " x 15", (11 cm x 38 cm)	1 ea
MDS91412LF	Infant	2 $\frac{7}{8}$ " x 10 $\frac{1}{4}$ ", (7 cm x 26 cm)	1 ea
MDS91413LF	Large Adult	6 $\frac{7}{8}$ " x 26", (17 cm x 66 cm)	1 ea
MDS91414LF	Thigh	8 $\frac{1}{2}$ " x 30 $\frac{1}{2}$ ", (22 cm x 77 cm)	1 ea

Two-Tube Seamless Inflation Bag

Item #	Description	Size	Pkg.
Latex			
MDS9505	Adult	4 $\frac{3}{4}$ " x 8 $\frac{7}{8}$ ", (12 cm x 23 cm)	1 ea
MDS9506	Large Adult	6" x 12 $\frac{3}{4}$ ", (15 cm x 32 cm)	1 ea
Neoprene Latex-Free			
MDS9505LF	Adult	4 $\frac{3}{4}$ " x 8 $\frac{7}{8}$ ", (12 cm x 23 cm)	1 ea
MDS9506LF	Large Adult	6" x 12 $\frac{3}{4}$ ", (15 cm x 32 cm)	1 ea

One-Tube Seamless Inflation Bag

Item #	Description	Size	Pkg.
Latex			
MDS9205	Adult	4 $\frac{3}{4}$ " x 8 $\frac{7}{8}$ ", (12 cm x 23 cm)	1 ea
MDS9206	Large Adult	6" x 12 $\frac{3}{4}$ ", (15 cm x 32 cm)	1 ea
Neoprene Latex-Free			
MDS9205LF	Adult	4 $\frac{3}{4}$ " x 8 $\frac{7}{8}$ ", (12 cm x 23 cm)	1 ea
MDS9206LF	Large Adult	6" x 12 $\frac{3}{4}$ ", (15 cm x 32 cm)	1 ea

Nylon Range finder Velcro Cuffs

Item #	Description	Size	Pkg.
MDS9420	Adult	5 $\frac{3}{4}$ " x 22", (15 cm x 56 cm)	1 ea
MDS9421	Child	4 $\frac{1}{4}$ " x 15", (11 cm x 38 cm)	1 ea
MDS9423	Large Adult	6 $\frac{7}{8}$ " x 26", (17 cm x 66 cm)	1 ea

ACCESSORIES

Gauges and gauge parts

Item #	Description	Pkg.
MDS9414	Premier Gauge Only	1 ea
MDS9415	Face Glass for Premier Gauge	1 ea
MDS9381	Standard Gauge Only	1 ea
MDS9507	Face Glass for Standard and Economy Gauge	1 ea



MDS9414



MDS9381

Bulbs and Valves

Item #	Description	Pkg.
MDS9508	Valve Only	1 ea
MDS91481	Bulb Only, Latex-Free (PVC)	1 ea
MDS91481LF	Bulb Only, Latex-Free (Neoprene)	1 ea
MDS91482	Bulb and Valve Assembly, Latex-Free (PVC)	1 ea
MDS9630	Bulb and Valve Assembly, Premium Latex	1 ea
MDS9560	Bulb and Valve, Large, Latex	1 ea



MDS9508



MDS9630



MDS91481LF

Tubing

Item #	Description	Pkg.
MDS9473	4-Foot Coiled Tubing w/Connector, Latex	1 ea
MDS9475	8-foot Coiled Tubing w/Connectors, Latex	1 ea
MDS9475LF	8-foot Coiled Tubing w/Connectors, Latex-Free (Neoprene)	1 ea



MDS9475LF

Other parts

Item #	Description	Pkg.
MDS9472	Inflation Systems Holder for Wall Mounted or Mobile Sphygmomanometers	1 ea
MDS9481	Metal Connector	1 ea



MDS9472



MDS9481

STETHOSCOPES

Our Best Value Stethoscopes!

Quality and economy makes our single and dual-head stethoscopes excellent values.

- Both stethoscopes feature 22" (56 cm) tubing and lightweight design.
- Thick-walled tubing minimizes background interference.
- Dual-heads feature a non-chill ring on the bell for increased patient comfort.

Item #	Description	Pkg.
Single-Head Stethoscopes		
MDS926101	Black	1 ea
MDS926102	Gray	1 ea
MDS926103	Blue	1 ea
MDS926105	Lavender	1 ea
MDS926106	Red	1 ea
MDS926107	Hunter Green	1 ea
MDS926108	Yellow	1 ea
Dual-Head Stethoscopes		
MDS926201	Black	1 ea
MDS926202	Gray	1 ea
MDS926203	Blue	1 ea
MDS926205	Lavender	1 ea
MDS926206	Red	1 ea
MDS926207	Hunter Green	1 ea



This Stethoscope has EXCELLENT Sound!

The Sprague Rappaport Stethoscope incorporates versatility into an excellent sound-transmitting device.

- Each stethoscope comes complete with a component kit that includes two sets of eartips, large and small diaphragm, adult bell, pediatric bell, and infant bell.
- The dual 22" (56 cm) tubing design transmits even the faintest sounds.

Item #	Description	Pkg.
MDS926301	Black	1 ea
MDS926302	Gray	1 ea
MDS926303	Blue	1 ea
MDS926305	Lavender	1 ea
MDS926306	Red	1 ea
MDS926307	Hunter Green	1 ea



Single or Dual-Head Parts

Item #	Description	Pkg.
MDS9544	White Plastic Ear Tips	1 ea
MDS9547	Soft Rubber Ear Tips	1 ea
MDS9552	Diaphragm	1 ea



Stethoscopes

STETHOSCOPES

Medline's Stainless Steel Stethoscopes offer unsurpassed quality and value

- Each stethoscope is handcrafted from high tensile strength German stainless steel. Features 22" (56 cm) tubing.
- The color coordinated non-chill rings provide excellent patient comfort.
- All models come in a foam-lined box with one extra set of ear tips.
- Available in adult, pediatric, neonatal, and cardiology styles.
- Cardiology Stethoscope has a deep, cone-shaped bell — a unique feature which provides excellent sound transmission and acoustic sensitivity. Features 17" (44 cm) tubing.

Item #	Description	Pkg.
Adult		
MDS92250	Gray	1 ea
MDS92260	Black	1 ea
MDS92270	Blue	1 ea
MDS92280	Burgundy	1 ea
Pediatric		
MDS92200	Gray	1 ea
Neonatal		
MDS92150	Gray	1 ea
MDS92160	Pink	1 ea
Cardiology		
MDS92500	Black	1 ea



Our Bowles Stethoscope features thick-walled 22" (56 cm) tubing, a chrome-plated binaural and a chrome-plated, raised, metal-stem chest piece.

Item #	Description	Color	Pkg.
MDS9515	Bowles Stethoscope	Black	1 ea

The Pediatric and Neonatal Stethoscopes have all the features of the dual-head stethoscopes, but are sized to fit the needs of the individual practitioner and patient.

- Both stethoscopes feature 22" (56 cm) tubing and lightweight design.
- Thick-walled tubing minimizes background interference.
- Aluminum binaurals.

Item #	Description	Color	Pkg.
MDS9557	Pediatrics	Gray	1 ea
MDS9565	Neonatal	Gray	1 ea



Stethoscopes

STETHOSCOPES



MIF767X04 – Navy Blue



MIF767XB0 – BlackOut

MDF® 2-in-1 Deluxe Sprague Rappaport Stethoscope

The advanced X-configuration acoustic tubing on the MDF® 767X eliminates all tube rubbing noises associated with the traditional Sprague-type stethoscopes and provides superior acoustic performance making it the most advanced Sprague Rappaport available today.

The chestpiece enables users to switch between adult and pediatric-sided bowels/bells with ease and interchange both sides of the chest piece with adult, pediatric and infant fittings (2 diaphragms and 3 bells). The 2-in-1 tubes in the X configuration converge at the middle of the tube and retain a dual lumen and dual input from the chestpiece. This innovation provides superior acoustic quality never available before on a Sprague Rappaport type stethoscope.

Item #	Color	Pkg.
MIF767X04	Navy Blue	1 ea
MIF767X08	Purple	1 ea
MIF767X10	Royal Blue	1 ea
MIF767X17	Burgundy	1 ea
MIF767XB0	BlackOut (All Black)	1 ea



MIF74010 – Royal Blue



MIF740B0 – BlackOut

MDF® Pulse Time™ Stethoscope

The original, patented MDF® 740 Pulse Time Stethoscope, integrated with a digital timing device, enables medical professionals to maintain eye contact with patients during diagnosis while monitoring elapsed time.

This innovative and patented Pulse Time stethoscope has a liquid-crystal analog (LCA) quartz timepiece equipped on top of the integrated bowels-type chestpiece.

The advanced chestpiece is made of lightweight anodized aluminum and the diaphragm is fitted with the unique, raised, ultra-thin fiber diaphragm for greater amplification and transmission of heart and lung sounds.

Item #	Color	Pkg.
MIF74008	Purple	1 ea
MIF74010	Royal Blue	1 ea
MIF74011	Black	1 ea
MIF74017	Burgundy	1 ea
MIF740B0	BlackOut (All Black)	1 ea

STETHOSCOPES

The Teaching/Training Stethoscope is designed for multi-function teaching applications.

- The dual head design offers the necessary flexibility to use this scope in all clinical and educational conditions.
- Comes complete with two binaurals, two pieces of Y-tubing, and an extension tube.

Item #	Description	Pkg.
MDS9750	Teaching/Training Stethoscope	1 ea



The Fetal Stethoscope is designed for superior sound sensitivity.

- Chrome-plated inner spring binaural.
- Vinyl coated headrest.
- 22" (56 cm) tubing.


Item #	Description	Pkg.
MDS9775	Fetal Stethoscope	1 ea



Our Disposable Stethoscope features plastic binaural and a sensitive chest piece that provides excellent sound quality.

- Designed for single patient use.
- Helps reduce cross-contamination in infectious areas.
- High visibility yellow.

Item #	Description	Pkg.
MDS9543	Disposable Stethoscope	1 ea

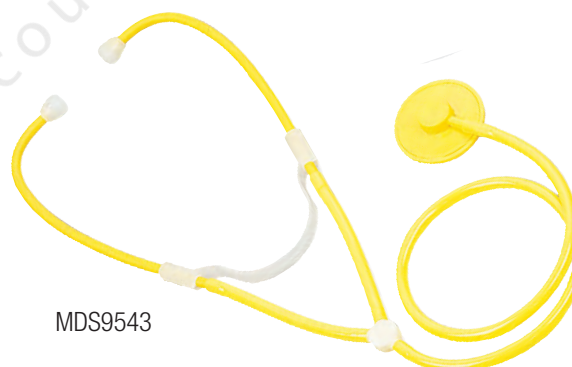

COMPASS DIABETES PROGRAM

Improving Diabetes Management in Long-Term Care

Medline's Compass Program A Diabetes Resource for Long-Term Care

Diabetes is a devastating disease affecting approximately 23.6 million people in the U.S. 5.7 million cases continue to go undiagnosed¹ for many years without treatment until severe complications occur. As many as 1 in 4 nursing home residents have diabetes.² The time and cost associated with caring for these individuals greatly impact the long term care facilities.

Medline has developed a comprehensive tool to help long-term care facilities provide better care for residents with diabetes by compiling and organizing data acquired during the course of working with healthcare professionals for more than 40 years. "Compass- A Diabetes Resource for Long-Term Care" emphasizes the importance of proper diabetes management to prevent complications and also illustrates how to achieve it.


Considering that heart disease, stroke, high blood pressure, blindness, kidney disease, dental disease, nervous system disease, and even amputations are all complications of diabetes¹, it is easy to see why healthcare providers strive to provide the best possible care for residents. Additionally, knowing that the Centers for Medicare & Medicaid Services is reviewing Pay for Performance measures and evaluating avoidable hospital admissions, like those due to complications from diabetes, it is vital for nursing home staff to understand current standards of practice.

With 26 percent of the nursing home population battling diabetes and its complications on a daily basis, Medline is providing educational tools for residents, their families and their caregivers.



References:

1. National Diabetes Fact Sheet, 2007. Centers for Disease Control and Prevention website. Available at www.cdc.gov/diabetes/pubs/pdf/ndfs_2007.pdf. Accessed on April 24, 2009
2. Johnson EL, Brosseau JD, Soule M, et al. Treatment of diabetes in long-term care facilities: a primary care approach. *Clinical Diabetes*. 2008;26(4):152-156

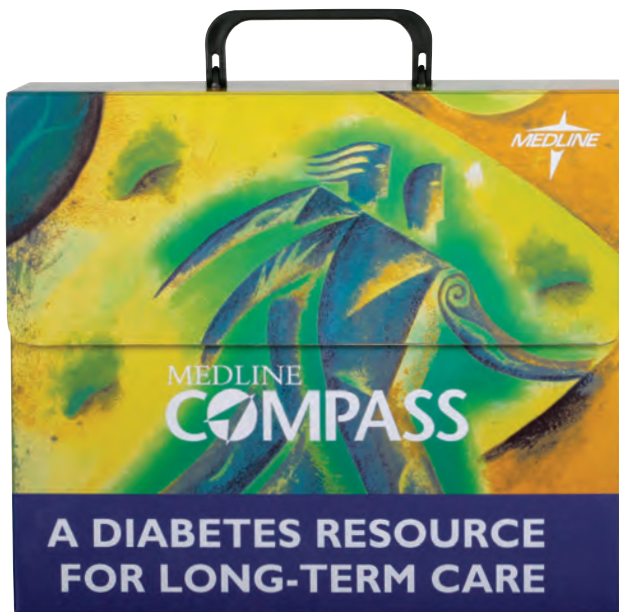


What's in the Compass Diabetes box?

- **One Director of Nursing (DON) Instruction Manual** — A teacher's guide to the program, which includes 4 CE contact hours.
- **Three additional Glucose Monitoring Training resource books** — Put them on treatment carts or in nurses' stations for easy access.
- **Five additional Family Education books** — Give them to residents' families.
- **One Educational Compact Disc and one Educational DVD**

To educate staff on using the OptiumEZ™ Glucose Monitoring System.

Compass Diabetes Program

A look at a few things from the program:

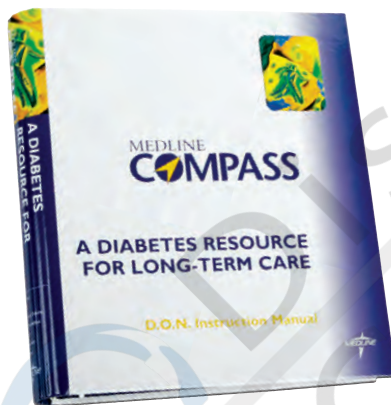


Instructional Manual



Educational CD and DVD

Resource Books

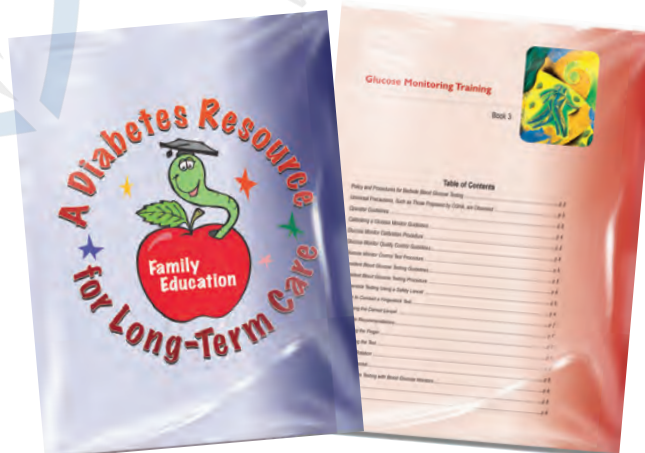

Family Educational Books

POINT OF CARE TESTING

New!



Point of Care Testing

New! CoaguSense PT/INR System

The CoaguSense PT/INR system provides a more direct approach to INR testing as the only portable monitoring system that directly detects clotting endpoint – providing greater confidence in results. Features convenient testing with no coding required, rapid and reliable results and a wide hematocrit range of 15 – 60%.

Item #	Description	Pkg.
ABO03P6001	CoaguSense PT/INR Professional System – Includes 1 Meter, 1 Carrying Case, 1 Professional User Manual, 1 Professional Quick Reference Guide, 1 Auto Single Usage Lancets – bag of 12, 4 AA alkaline batteries, 2 MiniPipettes (ABO03P5502Z), 1 Pipette Tips pack of 96 each (ABO03P5402Z)	1 ea
ABO03P5650Z	CoaguSense Test Strip Kit – 50 ea/bx, includes 54 transfer tubes (ABO03P5354Z). Must order 1 box of ABO03P7204Z at \$.01 with each box ordered	1 bx
ABO03P7204Z	CoaguSense Control Strip Kit – 4 ea/bx, includes 2 high and 2 low control strips and 1 vial of control solution. Must be ordered at \$.01 with each box of ABO03P5650Z	1 bx
ABO03P6910Z	CoaguSense Control Strip Kit – 10 ea/bx, includes 5 high and 5 low control strips and 1 vial of solution	1 bx
ABO03P5502Z	CoaguSense MiniPipette, pack of 2	1 pk
ABO03P5402Z	CoaguSense Pipette Tips, 96 each per tray	1 tray
ABO03P5354Z	CoaguSense Sample Transfer Tubes, pack of 54	1 pk

Professional Intended Use: The CoaguSense Prothrombin Time (PT)/INR Monitoring System is an in vitro diagnostic device that provides quantitative prothrombin time (PT) results, expressed in seconds and INR units. It uses fresh capillary whole blood. It is intended for use by health care professionals at the point of care to monitor patients who are on warfarin-type (Coumadin®) anticoagulation therapy. The device is not intended to be used for screening purposes.

POINT OF CARE TESTING

Medline OptiumEZ Glucose Monitoring System

The Medline OptiumEZ Blood Glucose Monitoring System requires no coding and provides easy, accurate, and reliable results!

Easy to Test

- No Coding Required
- Simple 2-step Testing
- Results in 5 seconds
- Small blood sample size – 0.6 uL

Easy, Accurate and Reliable Results

- Unique chemistry minimizes the potential interfering effects from many common substances*
- Allows for the ability to reapply more blood on the same test strip
- Individually wrapped test strips preserve the integrity of the test strips, which ensures an accurate glucose reading

Item #	Description	Pkg.
MIS71039Z	OptiumEZ Blood Glucose Meter	1 ea
MIS71041	OptiumEZ Glucose Test Strips	12 bx/cs, 50 ea/bx
MIS71041Z	OptiumEZ Glucose Test Strips	50 ea/bx
MIS71042	OptiumEZ Glucose Test Strips	6 bx/cs, 100 ea/bx
MIS71042Z	OptiumEZ Glucose Test Strips	100 ea/bx
MIS70550Z	MediSense Control Solution, 1 vial high, 1 vial low per box	1 bx

*Common interfering substances include acetaminophen, uric acid, ascorbic acid and salicylic acid. Do not use during xylose absorption testing. Study conducted in 2002 and 2005. Data on file, Abbott Diabetes Care Inc.



No Coding
Required!

ART22011-Rev-A.indd 1

Point of Care Testing

EvenCare® G2™ Blood Glucose Monitoring System

Advanced technology makes it easier to care for residents with diabetes.

- **Easy-to-use** – Voice assistance guides you through the testing process and provides you with your results.
- **Auto code system** – Minimizes errors and inaccurate results.
- **Fast test time** – Results in 6 seconds.
- **Small sample size** – Small blood requirements means less discomfort.
- **Alternate site testing** – Option to test on the forearm, palm or fingertips.
- **Memory option** – Stores 300 test times and can provide 7-, 14- and 30-day averaging.

Item #	Description	Pkg.
MPH1540	Glucose Meter	1 ea
MPH1550	Test Strips	50 ea/bx, 12 bx/cs
MPH1550Z	Test Strips	50 ea/bx
MPH1560	Control Solution High/Low	6 bx/cs
MPH1560Z	Control Solution High/Low	1 bx
MPHBAAAZ	Replacement Battery, AAA	24 ea/bx

New!



No Code

with Voice Guidance
"103 mg/dL"

LANCETS



DIB0310510

Hemoglobin A1C Test

The A1cNow® InView System offers incredibly fast, quantitative A1C results (in just 5 minutes) utilizing micro-optical technology that is integrated into a proprietary reusable monitor with single use cartridges. Accuracy is assessed by the same standards used by laboratories.

Item #	Description	Pkg.
DIB0310510	A1cNow InView Monitor	10/kit
DIB3021	A1c Now InView Monitor	20/kit

Lancets




MSC0002

Monofilament

This monofilament can be used as a routine test to assess a resident's foot sensitivity.

Item #	Description	Pkg.
MSC0002	Foot Monofilament, Diabetic Neuropathy Test	25/bx

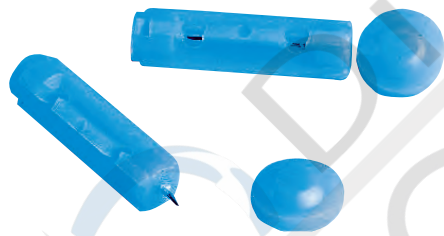

MPH100030

Lancets

Ultra Sharp Lancets have fine-point tip for smooth and fast procedures. These lancets can be used with all general lancing devices.

Item #	Description	Pkg.
MPH100030	Lancet, 30 Gauge	2500 ea/cs
MPH100030Z	Lancet, 30 Gauge	100 ea/bx



SEI100201Z

Lancing Device

Fits all general purpose lancets.

Item #	Description	Pkg.
SEI100201Z	Lancing Device, Adjustable, Non-Safety	1 ea

LANCETS

Medline Safety Lancets


Medline safety lancets are simple and safe for both the healthcare professional and the patient. The needle retracts automatically and locks out for extra safety, protecting patients and healthcare professionals from needlestick injuries and reducing cross-contamination.

- Needle retracts automatically and locks out for extra safety, protecting patient and health professionals from needle stick injuries.
- Different sizes available allowing health professionals to derive the proper amount of blood for the capillary test.
- Quick and simple to use, requiring little instruction.
- Cannot be reloaded or reused, protecting both patients and health professionals.
- Double spring action to provide greater patient comfort by firing and retracting the lancet with spring action.
- Designed to provide controlled needle depth allowing for greater patient comfort.
- Ergonomic and easy to handle with side ribs allowing for safe handling.

Item #	Description	Pkg.
MPHSAFETY21	Lancet, Normal Flow, 21 Gauge 1.8 mm	18 bx/cs, 200 ea/bx
MPHSAFETY21Z	Lancet, Normal Flow, 21 Gauge 1.8 mm	200 ea/bx
MPHSAFETY28	Lancet, Micro Flow, 28 Gauge 1.6 mm	18 bx/cs, 200 ea/bx
MPHSAFETY28Z	Lancet, Micro Flow, Gauge 1.6 mm	200 ea/bx



MPHSAFETY28

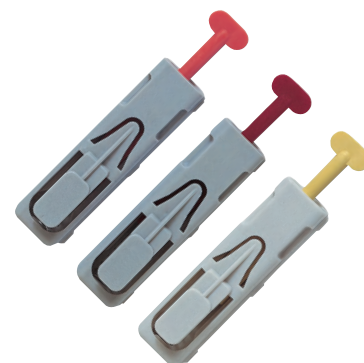

OMMAT07 series



OMMAT1004Z

Unistik® Safety Lancets

Offers maximum comfort for your patients. The Unistik needlepoint is hidden before use and automatically retracts after use to avoid needle stick injuries and help reduce cross-contamination.

Item #	Description	Pkg.
OMMAT0747Z	Unistik 2 Gentle	200/bx
OMMAT0704Z	Unistik 2 Normal Depth	200/bx
OMMAT0704	Unistik 2 Normal Depth	6000/cs
OMMAT0714Z	Unistik 2 Extra Depth	200/bx
OMMAT1004Z	Unistik 3 Normal Depth	200/bx

SHARPS CONTAINERS



MDS705153 Shown with
SWD85161H cabinet



MDS705151



MDS705203



MDS705115



MDS705110 shown
with Bracket MDS707912

Medline has a variety of sharps disposal systems to fit your needs.

And if you're looking to cut costs, you'll find savings, too, without having to sacrifice quality. Medline Sharps Containers are among the safest containers on the market. Easy-to-see fill levels reduce overfilling and accidental needle sticks. Containers lock easily for final disposal and safe transportation. Specially designed lids limit unwanted access to container contents.

Patient Room Sharps Containers

These sharps disposal containers are designed for use in patient rooms, on med carts and in high traffic areas.

- Lid design facilitates one-handed operation
- Horizontal drop design for maximum capacity to fill line
- Easy-to-see fill levels
- Large openings accommodate larger syringes
- Multiple cabinet and bracket options

Item #	Description	Pkg.
MDS705125	2.5 Quart (2.4 L), Wall Mount 7.5" H x 4.25" L x 11" W (191 x 108 x 279 mm)	30/cs
MDS705153	5 Quart (4.7 L), Wall Mount 11" H x 4.75" L x 10.75" W (279 x 121 x 273 mm)	20/cs
PLA141020	5 Quart (4.7 L), Wall Mount / Counter- Balanced Lid, 10.75" H x 4.75" L x 10.5" W (273 x 121 x 267 mm)	20/cs
MDS705151	5 Quart (4.7 L), Wall Mount / Counter- Balanced Lid, 14" H x 6" L x 13" W (356 x 152 x 330 mm)	14/cs
MDS705203	3 Gallon (11.4 L), Wall Mount/ Counter- Balanced Lid, 20.5" H x 6" L x 14" W (521 x 152 x 356 mm)	12/cs

Phlebotomy Sharps Containers

Portable sharps containers are designed to fit into blood drawing trays for point-of-use disposal. Square base for added stability.

- Translucent top allows easy identification of contents
- Choice of single or dual openings to accommodate a variety of sharps
- Containers lock for final disposal and safe transport
- Brackets and holders available

Item #	Description	Pkg.
MDS705110	1 Quart Phlebotomy Container (0.9 L), 6.25" H x 4.25" L x 4.5" W (159 x 108 x 114 mm)	100 ea/cs
MDS705115	1.5 Quart (1.4 L) Transportable Container 10.5" H x 3.5" L x 3.5" W (267 x 89 x 89 mm)	20 ea/cs

SHARPS CONTAINERS



MDS705210



MDS705201

Multipurpose Sharps Containers

These containers are designed for use in restricted access areas.

- Containers lock for final disposal
- Nestable containers save space and reduce shipping costs
- Non-tortuous lid design accommodates a variety of sharps sizes
- Horizontal drop maximizes use of container volume
- Brackets available

Sharps / Biohazard

Item #	Description	Pkg.
MDS705201	1 Gallon (3.8 L) 6.25" H x 7" L x 10.5" W (159 x 178 x 267 mm)	32/cs
MDS705202	2 Gallon (7.6 L) 10.25" H x 7" L x 10.5" W (260 x 178 x 267 mm)	20/cs
MDS705210	10 Gallon (37.9 L) 21.5" H x 12" L x 15.5" W (546 x 305 x 394 mm)	6/cs

Mail-Back Sharps Containers

Each kit includes all components required for collection and transportation of sharps and regulated medical waste.

- Compact design uses a minimum of counter space
- Pre-paid, pre-addressed shipping labels and pre-paid disposal
- Lids provide temporary and final locking closure
- Easy-to-use enclosed outer shipping box
- Step-by-step instructions
- Tracking documents

Item #	Description	Pkg.
SCC1QV2	1 Quart (0.9 L)	1 ea
MDS705000	1 Gallon (3.8 L)	1 ea
SCC2GV2	2 Gallon (7.6 L)	1 ea
SCC3GV2	3 Gallon (11.4 L)	1 ea
SCC5GV2	5 Gallon (18.9 L)	1 ea
SCC18GV2	18 Gallon (68.1 L)	1 ea



SHARPS CONTAINERS

Sharps Containers and Accessories at-a-Glance

Item #	Description	Locking Bracket	Non-Locking Bracket	Cabinet	Glove Box for Cabinet	Holder	Cart
MDS705110	1 Quart Phlebotomy, Red		MDS707912			MDS707910	
MDS705115	1.5 Quart Transportable, Red						
MDS705153	5 Quart, Wall Mount, Red	SWD8518X	SWD8519C	SWD85165H or SWD85161H	SWD8555SA		
PLA141020	5 Quart (4.7 L), Wall Mount/Counter-Balanced Lid, Red			SWD85165H or SWD85161H	SWD8555SA		
MDS705151	5 Quart, Wall Mount, Counter-balanced Lid, Red	MDS707054	MDS707062	MDS707553	SWD307013 or SWD320032	MDS707955	
MDS705203	3 Gallon, Wall Mount, Counter-balanced Lid, Red	MDS707054	MDS707062	MDS707005	SWD307013 or SWD320032	MDS707016	
MDS705201	1 Gallon, Sliding Lid, Red	MDS707963	MDS707975			MDS707974	
MDS705202	2 Gallon, Sliding Lid, Red	MDS707963	MDS707975			MDS707974	
MDS705210	10 Gallon, Hinged Lid, Red						MDS707109

Medication Disposal System

This is a complete disposal solution for unused medications. The medication-disposal system includes the cost of the medication-disposal container, return freight and total destruction costs. The unused medications can be over-the-counter (OTCs) drugs or prescription drugs, so long as the prescription drugs are not controlled substances.

Item #	Description	Pkg.
MDS707001	Medication Disposal, Envelope 8" x 11"	1 ea
MDS707010	Medication Disposal, Container 10 Gallon (37.9 L)	1 ea
MDS707020	Medication Disposal, Container 20 Gallon (75.7 L)	1 ea





SWD511110Z

SHARPS CONTAINERS

Safety Insulin Syringes

These safety syringes have a secure safety shield that covers the entire needle length.

Item #	Description	Pkg.
SWD511110Z	1 mL syringe, 29 gauge x 1/2" needle	100/bx
SWD511110	1 mL syringe, 29 gauge x 1/2" needle	500/cs
SWD511136Z	1/2 mL syringe, 29 gauge x 1/2" needle	100/bx
SWD511136	1/2 mL syringe, 29 gauge x 1/2" needle	500/cs
SWD511144Z	3/10 mL syringe, 29 gauge x 1/2" needle	100/bx
SWD511144	3/10 mL syringe, 29 gauge x 1/2" needle	500/cs

Medline's Tympanic Thermometer with Quick Probe Release

Get fast, accurate and easy-to-read results with Medline's Tympanic Thermometer with Quick Probe Release. This ear thermometer is designed as a comfortable and convenient method for taking a person's body temperature. No embarrassing or invasive procedures; you simply place the thermometer in the ear to get a reading within seconds. It's perfectly suited for both home and clinical use on patients of all ages. Plus, its unique mechanical probe release helps reduce cross-contamination, allowing you to easily change probe covers with every use!

Why an ear thermometer?

It turns out that the eardrum is an extremely accurate point to measure body temperature from because it is recessed inside the head (just like your tongue). Instead of touching the eardrum with a thermometer, the tympanic thermometer uses an infrared sensor to measure the temperature of energy radiating from the eardrum.

Innovative, quick-response ear thermometer features:

- One-touch operation
- Mechanical probe cover quick-release allows easy probe cover change
- Includes 100 probe covers
- Safe to use on patients of any age
- Readings in Fahrenheit and Celsius
- Takes an accurate, non-invasive temperature from the tympanic membrane

Item #	Description	Pkg.
MDS9700	Thermometer	1 ea
MDS9701	Probe covers	100/bx



Apply disposable cover onto the tympanic thermometer.

Turn device "on" by pressing button on front.



Place thermometer in ear and depress top button to take reading.

After reading, flip the blue tab to release the disposable probe cover without contact.



PATIENT MONITORING

Vital Signs Monitoring Solutions

247™ Diagnostic Station An Expandable Diagnostic Station with BP, Temperature and SpO₂

Medline is always looking for new and innovative products that bring unique value to the healthcare industry. We are excited to introduce the 247 Diagnostic Station from SunTech Medical®. SunTech has a 20-year history of research in blood pressure technology, and their pragmatic approach to blood pressure measurement was a good match with the common-sense approach we bring to our customers.

The new SunTech 247 Diagnostic Station combines the efficiency of an automated monitor with the convenience of a manual mode that allows the clinician to take a BP reading with a stethoscope whenever needed. Also, the 247 diagnostic station is fully modular; you can start with BP monitoring only, then add SpO₂ monitoring, or patient temperature or both.

In healthcare settings, the time-saving efficiency of an automated reading is highly desirable. However, there are patient conditions (arrhythmia), patient groups (pediatrics), and environments (moving vehicles) where automated monitors do not always function well.³ In these situations, having a manual mode built into the automated device is eminently practical.



Item	Description	Pkg.
Mobile Stand		
SNT99006200	BP	1 ea
SNT99006300	BP & Temperature	1 ea
SNT99006400	BP & SpO ₂	1 ea
SNT99006100	BP, Temperature & SpO ₂	1 ea
Wall Mount System		
SNT99006600	BP	1 ea
SNT99006700	BP & Temperature	1 ea
SNT99006800	BP & SpO ₂	1 ea
SNT99006500	BP, Temperature & SpO ₂	1 ea
Table Top System with Battery		
SNT99007000	BP	1 ea
SNT99007100	BP & Temperature	1 ea
SNT99007200	BP & SpO ₂	1 ea
SNT99006900	BP, Temperature & SpO ₂	1 ea
Modules		
SNT98012800	Temperature Module	1 ea
SNT98012900	SpO ₂ Module	1 ea

Accessories Included

BP: 8' (2.5 meters) blood pressure hose, Adult & Large Adult all purpose cuffs, mains/power supply, mains/power cable, resource CD, quick start guide, reference card

BP + Temperature: oral/axillary probe (blue), 1 box of Filac® FasTemp probe covers

BP + SpO₂: reusable adult finger sensor, 6' (1.8 meters) extension cable

Additional Accessories (not included)

Item	Description	Pkg.
SWD202020Z	Probe Covers for the FastTemp Thermometer	
STCUS1320MQ	Child Replacement Cuff	5/bx
STCUS2635MQ	Adult Replacement Cuff	5/bx
STCUS3242MQ	Large Adult Replacement Cuff	5/bx
STCUS4250MQ	Thigh Replacement Cuff	5/bx

Reference:

3. Jones DW, Appel LJ, Sheps SG, Rocella EJ, Lenfant C. Measuring blood pressure accurately: new and persistent challenges. *Journal of the American Medical Association*. 2003; 289(8):1027-1030.

CASMED 740 Vital Signs Monitor

Changes in the hospital environment are driving the need for more flexible solutions. The CASMED portable vital signs monitor is capable of spot-check and continuous monitoring, regardless of where the patient is located or how often they are transported. The monitor is rugged, portable and lightweight (weighs only 3 pounds). It is reliable and easy-to-use with multi-parameter capabilities, including:

- Fast, accurate blood pressure readings for neonates, children and adults. Motion does not affect readings with the CASMED.
- Monitoring of oximetry and temperature continuously or in "spot-check" mode

The CASMED 740 Vital Signs Monitor is also ideal for medical-surgical patients on pain medications (including PCA pumps), for hospital or office-based endoscopy procedures, and for use onboard ambulances and paramedic vehicles.



CASMED® 740



STC01020800

Item #	Description	Pkg.
STC01020800	NIBP Monitor	1 ea
STC01020801	NIBP & Masimo SpO2	1 ea
STC01020802	NIBP & Nellcor SpO2	1 ea
STC01020803	NIBP & Temp	1 ea
STC01020804	NIBP, Masimo SpO2 & Temp	1 ea
STC01020805	NIBP, Nellcor SpO2 & Temp	1 ea

Mounting Solutions

STC01020172	Stand with basket	1 ea
STC01020243	Universal mount 2 hole	1 ea
STC01020245	Universal mount 4 hole	1 ea
STC01020173	Swivel mount	1 ea
STC01020181	Printer bracket	1 ea
STC01020265	12" wall mount	1 ea
STC01020447	Mounting adapter	1 ea
STC01020448	Mounting plate	1 ea
STC01020189	External Printer	1 ea

THERMOMETERS



TMS40002

Thermoscan® Thermometers

The Braun® Thermoscan Pro 4000 Thermometer patented technology is an ergonomically designed ear thermometer engineered to take temperatures easily and accurately with the click of a button.

- Exclusive ExacTemp® technology provides an active user feedback light that signals that the probe has been positioned in a stable way through out the measurement, avoiding inaccurate measurements by handling errors
- Preheated probe tip minimizes the cool down effect experienced when taking repeat temperatures, improving accuracy and repeatability
- Automatic probe cover eject button for quick and easy removal of used probe covers
- Built-in cross-contamination safeguards – probe cover sensor reminds operator to use a new probe cover
- Wide range of patient temperature – 68°F-108°F
- Wide range of ambient temperature – 50°F-104°F ambient range allows for use in most environments
- Durable – tested for damage resistance
- Convenient wall-holder for Pro 4000 unit w/probe cover storage is available
- Three-year limited warranty coverage

Item #	Description	Pkg.
TMS40002	Thermoscan PRO 4000	1 ea
TMSPC800	Probe Covers	40 bx/cs



MLI111210

NexTemp® Disposable Clinical Thermometers

- Clinical Accuracy
- Helps to Prevent
- Cross-Infection
- Resets Automatically
- No Special Storage
- Requirements
- Carries FDA, ISO & CE Certifications
- 5-Year Shelf Life
- Latex-Free
- Non-Toxic
- Mercury-Free
- Economical

Item #	Pkg.
NexTemp Fahrenheit	
MLI111210	100/bx
MLI111220	20 bx/cs
NexTemp Celsius	
MLI211210	100/bx
MLI211220	20 bx/cs

THERMOMETERS

Exergen® Temporal Scanner™ Thermometer

Precision infrared scanner provides a truly non-invasive method for measuring temperature. Simply scan across the temporal artery area of the forehead, touch behind the earlobe, and an accurate reproducible temperature is measured in about three seconds.

- Convenient one-handed operation
- Completely non-invasive
- Stainless steel probe for highest durability
- Easy snap-on probe caps
- Three-year warranty

Item #	Description	Pkg.
EXGTAT5000	Infrared Temporal Artery Thermometer	1 ea
EXGTAT2000	Infrared Temporal Artery Thermometer	1 ea
EXG134203	Disposable Probe Caps	1000 ea/cs



EXGTAT5000

Six Seconds*, Temple Digital Thermometer

The Medline Instant Read* Digital Temple Thermometer accurately measures body heat using Patented R.A.T.E.™ technology (Rapid, Accurate, Temperature, Establishment).

- Sensors rapidly track heat flow generated from blood vessels to the skin's surface and convert measurements to body temperature.
- Easy-to-use while asleep or awake.
- Meets accuracy standards.
- Large digital display.
- Switchable from Fahrenheit to Celsius.
- Beeps when ready.
- Memory stores previous reading.

Item #	Description	Pkg.
MDS9698	Temple Thermometer	1 ea
MDS9699	Probe Covers	250 ea/bx

*Readings in six to ten seconds depending on individual.



Gentle, non-invasive sensor accurately reads in F° & C°

Power/memory control button

Large Digital Display

2-AAA Batteries included

Features Patented R.A.T.E. Technology



Thermometers

THERMOMETERS



SWD504000

Get Both Speed And Accuracy with the FILAC® 3000 Features Both!

The FILAC 3000 Electronic Thermometer utilizes advanced technology to ensure the most accurate and fastest readings are provided in easy-to-use formats. These cutting edge thermometers address the need for fast, accurate readings in under 4 seconds or an advanced accuracy reading in under 6 seconds. Compact, user-friendly design will stand upright on a flat surface; intuitive icon-based interface system with backlight LCD display; unique progress bar indicates status of temperature measurement; last temperature recall function. Patented isolation chamber design and color schemes along with the durable probe covers help minimize the chance of cross-contamination and enhance infection control measures. Operates on 4 AA batteries for 6000 measurements. Equipped with electronic anti-theft system, self-calibrating electronics and auto on-off function. Includes oral probe and push button probe cover ejection.

Item #	Description	Pkg.
SWD504000	Filac 3000 Thermometer	1 ea
SWD202020	Probe Covers for the Filac Fastemp Thermometer & 3000	5,000/cs
SWD202020Z	Probe Covers for the Filac Fastemp Thermometer & 3000	500/box



SWD303000

Genius2 Thermometer

Measures temperature from the tympanic membrane, which shares the same blood flow as the hypothalamus, the body's natural thermostat.

- Utilizes infrared technology to derive the most accurate body temperature from the tympanic membrane.
- Records and analyzes up to 100 temperature measurements, and selects the highest and most accurate in less than 2 seconds, as compared to 32 in Genius and 1-8 measurements in competitive thermometers.
- Reduces the chance of an inaccurate reading due to improper technique or patient movement.
- Durable tip-based sensor technology Improves durability and accuracy, minimizing the need for recalibration.
- Tapered probe design
- Rapid response time 1 - 2 second
- Real time temperatures
- Programmable display
- Reduces cross-contamination
- No contact with mucus membranes
- No-touch probe cover application and release

Item #	Description	Pkg.
SWD303000	Genius2 Thermometer	1 ea
SWD303030Z	Probe Covers	96 ea/bx
SWD303030	Probe Covers	2,112/cs

THERMOMETERS

Quick. Accurate. Reliable

Infrared Thermometer with Digital Readout

- Professional accuracy
- Fast 1-second readout
- Wide range readout 50°F – 122°F person 32°F – 221°F object/liquid
- Nite-Glo® display
- Fever alarm
- 10-reading memory recall
- °F or °C switchable
- Automatic shut-off

- Includes:**
- Convenient storage bag
 - Detailed guidebook in English and Spanish
 - 2 AAA batteries
 - 5-year limited warranty



Dual Scan Modes

Forehead and object/liquid

Easy-To-Read

Nite-Glo® display

Item	Description	Pkg.
MBH18535000	Infrared Thermometer	1 ea

MBH18535000

DeNovo Broadband Liquid Crystal Temperature Indicator

This disposable forehead thermometer provides the widest temperature range at the lowest price: 90 – 106 degrees F and 33– 41 degrees C.

- Clearly marked “normal” core temperature
- Measures in two-degree increments
- Latex-free



DVPBB106

Item #	Description	Pkg.
DVPBB106	Disposable forehead thermometer, 2 degree increments, 90°F – 106°F and 33°C – 41°C	100/bx

DeNovo Temp Assure Temperature Trend Indicator

This disposable forehead thermometer provides temperature ranges from 94 - 105 degrees F and 34– 40 degrees C.

- Clearly marked “normal” core temperature
- Measures in one-degree increments
- Allows for patient monitoring during transport
- Great for use with laryngeal mask airways
- Latex-free



DVPBB105

Item #	Description	Pkg.
DVPTA105	Disposable forehead thermometer, single degree increments, 94°F – 105°F and 34°C – 40°C	100/bx

Thermometers

THERMOMETERS




Digital Thermometers

Our digital thermometers provide accurate temperatures in about 60 sec.

- Large LCD screen for easy readings.
- Simple push-button start.
- Disposable sheaths allow for sanitary oral, rectal, or axillary use.

Item #	Description	Pkg.
MDS9650	Premier Oral °F & °C	1 ea
MDS9652	Premier Rectal °F	1 ea
MDS9654B	Standard Oral °F	200/cs
MDS9628	Oral Kit w/20 Sheaths °F	1 ea
MDS9629	Rectal Kit w/20 Sheaths °F	1 ea
MDS9607	Oral Sheaths	100/bx
MDS9608	Rectal Sheaths	100/bx
MDS96801	Flex-tip, waterproof, °C or °F	1 ea
MDS9655	Standard, Oral, Bulk, °F & °C	200/cs
MDS9655Z	Standard, Oral, °F & °C	20/bx
MDS9650B	Premier Oral, Bulk, °F & °C	200/cs
MDS9656	Standard, Rectal, Bulk, °F & °C	200/cs
MDS9656Z	Standard, Rectal, °F & °C	20/bx

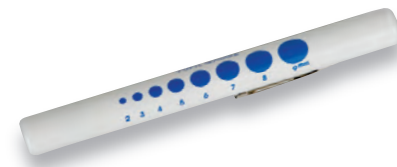

Disposable Penlights

Lightweight and disposable, these high-intensity exam lights are perfect for the physician, nurse, and EMT. They're constructed of high-impact molded styrene and can be clipped on pockets.

Item #	Description	Pkg.
MDS131040	Disposable Penlight, White	6/pk

DEFIBRILLATORS

Be Prepared with the Zoll AED PLUS® External Defibrillator

Providing high-quality cardiopulmonary resuscitation (CPR) and defibrillation are important components in helping to save someone's life from sudden cardiac arrest. Only ZOLL® offers innovative Real CPR Help® technology, providing real-time feedback for both rate and depth of chest compressions. The advantage of the AED Plus is that it may be used by anyone who has received basic training on its operation. 7 Year warranty on the device.

Item #	Description	Pkg.
ZOL21110201PK	ZOLL AED PLUS® Kit Includes: 1 Semi Automatic AED Plus Unit, 1 wall Cabinet, 1 Set of CPR-D PADZ, 1 Set of STAT II PADZ, 1 Set of Batteries (Also available in carrying case option (no wall mount) – item # ZOL21110201CC)	1 ea
ZOL8900080001	CPR-D PADZ ADULT replacement pads One piece defibrillation and CPR System with compression, depth and rate sensors. Supplied with gloves, barrier mask, scissors, razor, wet wipe and dry wipe. Five (5) year shelf-life.	1 ea
ZOL89000801	ADULT STAT PADZ II replacement pads 1 pair of multifunction electrodes	1 ea
ZOL8000080701	Lithium Battery Type 123 Lithium Batteries, quantity of 10 with storage sleeve	1 ea



ZOL21110201PK

Powerheart® AED G3 External Defibrillator

Cardiac Science® PowerHeart AED G3 Automated External Defibrillator with voice prompts and text screen displays guide you through the rescue process. When cardiac arrest occurs, every second counts. The advantage of the PowerHeart AED G3 is that it may be used by anyone who has received basic training on its operation. Electrodes are pre-connected and disposable. Semi-automatic defibrillator with biphasic variable waveform delivers an effective shock in less than 55 seconds. Daily, weekly, and monthly automated self-checks ensure that PowerHeart is always ready to go. Internal memory holds 60 minutes of ECG data. Includes one pair of defibrillation electrodes, instructional CD-ROM and training video, RescueLink® and MDLink® software and serial communication cable. Seven-year warranty on device & four-year full replacement on battery.

Item #	Description	Pkg.
AED9390SEMICC	Semi Automatic AED Kit Includes: 1 Replacement set of adult defibrillator electrode pads, 1 Soft carrying case, 1 Ready Kit. Also available with wall mount (item # AED9390SEMIPK)	1 ea
AED9390AUTOCC	Fully automatic AED Kit Includes: 1 Replacement set of adult defibrillator electrode pads, 1 Soft carrying case, 1 Ready Kit. Also available with wall mount option (item # AED9390AUTOPK)	1 ea



Powerheart G3 AED

Defibrillators

BioCon-500 Scanner safely measures bladder volume.

The BioCon-500 Scanner is brought to you exclusively by Medline

In our commitment to bring you high quality products at a good value, Medline has become the exclusive distributor of the BioCon-500—the only portable, 3-D ultrasound device with a pre-scan feature. The BioCon-500 is a quick, safe and non-invasive way to measure bladder volume and post-void urinary retention.

Help save time, increase accuracy with our bladder scanner

Incontinence is widespread among the elderly population, especially those who live in nursing homes, where 50 percent or more of residents are incontinent.⁴ The BioCon-500 can be helpful in developing individualized urinary continence programs.

As part of a urinary continence program, the BioCon-500 Bladder Scanner allows caregivers to:

- Scan the resident's bladder to determine if catheterization or toileting is necessary.
- Show residents the results so they can actually see how much fluid is in their bladder. This helps residents re-learn the sensation of needing to urinate.
- Develop individualized toileting schedules based on each resident's bladder function pattern rather than at fixed times, such as every two hours. This will save staff valuable time.
- More accurately determine if a resident's bladder volume and voided amounts are consistently low, prompting an increase in the resident's fluid intake.

Not for fetal use or for use on patients who are pregnant.



Item #	Description	Pkg.
MDSB10500	BioCon-500 Bladder Scanner	ea
MDS3110100032	Mobile Stand	ea
MDS092005	Ultrasound Gel	12 ea/cs

4. U.S. Department of Health & Human Services. Agency for Healthcare Research and Quality. Urinary Incontinence in Adults: Clinical Practice Guideline Update. Available at: <http://www.ahrq.gov/clinic/uiovervw.htm>. Accessed September 22, 2009.

ONLINE CATALOG NOW AVAILABLE

- 1 Download a QR Code Reader app
- 2 Launch the QR app
- 3 Scan this QR Code or visit <http://www.medline.com>



Medline Industries, Inc.

One Medline Place
Mundelein, IL 60060

Medline United States
1-800-MEDLINE (633-5463)

www.medline.com | info@medline.com

Medline Canada
1-800-396-6996

www.medline.ca | canada@medline.com

Medline México
01-800-831-0898

www.medlinemexico.com | mexico@medline.com

©2011 Medline Industries, Inc. A1cNow is a registered trademark of Metrika, Inc. Braun is a registered trademark of Braun GmbH Corporation. Cardiac Science and Powerheart are registered trademarks of Cardiac Science Corporation. Coagusense is a trademark of CoaguSense, Inc. -ExacTemp is a registered trademark of Braun GmbH Corporation. Exergen is a registered trademark and Temporal Scanner is a trademark of Exergen Corp. Fastemp is a trademark of Westek Corporation. Feverscan is a trademark of LCR Hallcrest. FILAC is a registered trademark of Covidien AG Corp. INRatio is a registered trademark of HemoSense, Inc. MDLink is a registered trademark of Cardiac Science Corporation. MDF is a registered trademark and PulseTime is a trademark of MDF Instruments USA, Inc. Medline is a registered trademark of Medline Industries, Inc. Nite-Glo is a registered trademark of MABIS Healthcare, Inc. OptiumEZ is a trademark of Abbott Laboratories. PowerHeart is a registered trademark of Cardiac Science Corp. R.A.T.E. is a trademark of BestMed. RescueLink is a registered trademark of Cardiac Science Corporation. Suntech is a registered trademark of Suntech Medical. Thermoscan is a registered trademark of Braun GmbH Corp. Tyvek is a registered trademark of E. I. Du Pont De Nemours & Co. Unistik is a registered trademark of Owen Mumford Ltd. Welch Allyn and SureTemp are registered trademarks of Welch Allyn, Inc. ZOLL and AED Plus are registered trademarks of ZOLL Medical Corporation. Medline, EvenCare, and NexTemp are registered trademarks and G2 and Simple Solutions. Significant Results are trademarks of Medline Industries, Inc.

MKT211181 / LIT007A / 5M / SG25