



SONICAID ENCORE

FETAL / MATERNAL INTRAPARTUM MONITOR



Outstanding Performance

The high performance Encore fetal and fetal / maternal monitors use today's most advanced technology to accurately acquire data for clinician decision support. Whether you need standard fetal monitoring, or more advanced fetal and maternal care, the Encore monitors provide accuracy and clarity when you need them.



Advanced Technology

"Widebeam" ultrasound technology provides a wide, deep, uniform beam shape which, coupled with high sensitivity and superb digital Doppler audio, enables clinicians to quickly and easily locate and maintain contact with the fetal heartbeat.

Encore incorporates unique technology to ensure accurate tracking in the most demanding of circumstances, minimizing the requirement for direct fetal ECG monitoring:

- 'Locate and Track' dynamic range gating technology locks onto the fetal heart signal for optimum pick-up.
- Advanced fetal heart rate algorithms with artifact rejection for outstanding FHR performance.

Excellent low-perfusion **SpO₂ performance**. Maternal heart rate displayed and printed as a trace for quick and easy differentiation from the fetal heart rate trace.

Quiet and fast **NIBP operation** minimizes distraction and discomfort for mom.



High Sensitivity Toco Transducer

The advanced design and high sensitivity of the toco transducer enable uterine activity monitoring for patients of all sizes, including those with high BMI.





Attention to Detail


- Slimline profile for maximum patient comfort
- Waterproof for easy cleaning (not for underwater use)
- Color coded, matching connectors (plugs and sockets) for intuitive use
- Built-in rack provides storage for transducers not in use





Serial and LAN interfaces provide seamless integration into most central recording and OB data management systems.

Color display makes it easy to monitor both maternal and fetal vital signs at a glance. Transducers match on-screen colors for fast, error free recognition.

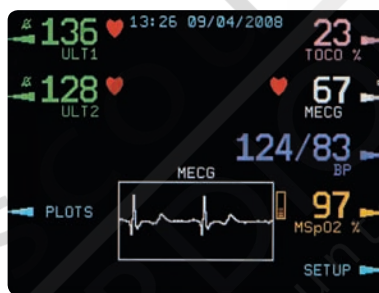

Capability Summary	94820	94830
Single or Twins Ultrasound FHR	✓	✓
Highly Sensitive Toco	✓	✓
Internal Uterine Activity	✓	✓
Direct Fetal ECG	✓	✓
Maternal ECG	✓	✓
Maternal SpO ₂		✓
Maternal NIBP		✓

- Ergonomic design fits into most cabinets, creating a soft and unobtrusive presence.
- Optional slim-line roll stand is easy to maneuver with integrated supplies tray.
- Wall mount also available.

www.spacelabshealthcare.com



Spacelabs Healthcare Company Headquarters
5150 220th Ave. SE, Issaquah, WA 98029
PO Box 7018, Issaquah, WA 98027-7018

Telephone: +1 425-657-7200
In North America: +1 800-522-7025
Fax: +1 425-657-7212

© 2011 Spacelabs Healthcare
Specifications subject to change without notice. All brands and product names are trademarks of their respective owners.

Product may not be available in all locations. Check with your local Spacelabs representative.

M1122-my11