



Barrier-Free® Power Examination Table

625





## Easier access puts patients and caregivers in their comfort zone



The Midmark 625 Barrier-Free® Examination Table – the most patient-accessible table available today.

Whether transferring a wheelchair patient or assisting an expectant mom, one thing is clear: the Midmark 625 means that positioning patients onto the exam table will be easier, more efficient, and much safer for you and your staff. With its 18- to 37-inch adjustable table height, lifting is less stressful, more controlled, and helps ensure a patient's dignity in the exam room.



## Capacity, connectivity, comfort... and more

As healthcare continues to evolve, it's essential to keep pace so providers and caregivers can move confidently into the future. There are enormous benefits to be gained by streamlining the exam room workflow, and one reliable way to achieve this is with the Midmark 625 Barrier-Free® Examination Table. The Midmark 625 offers all the capacity, accessibility and connectivity needed to provide excellent patient care without additional strain on you or your staff. Plus, with timesaving new features built into the Midmark 625, you'll be able to face the future of healthcare knowing that time truly is on your side.

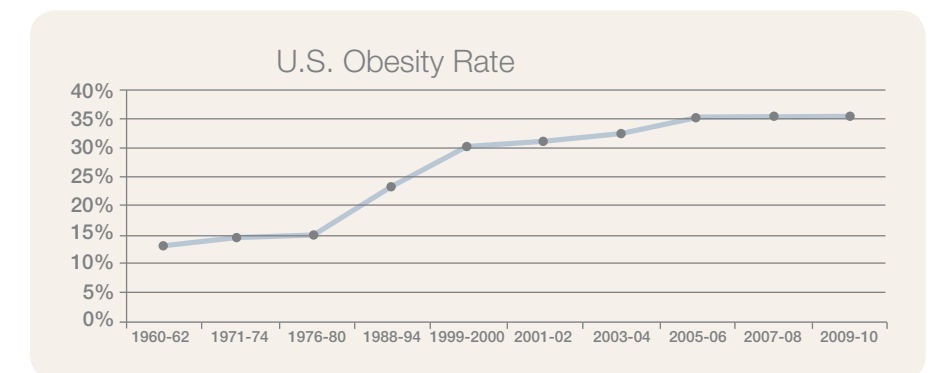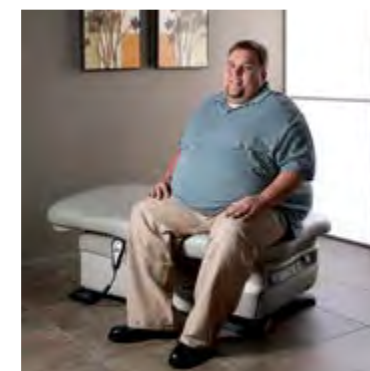

Source: CDC/NHAINES: adults ages 20-74 years

### Support for a changing demographic

Over one third of the population is considered obese, and the number continues to rise. Many healthcare providers face the risk of being ill equipped to care for a growing number of patients.

Don't let your equipment put limitations on your practice. The Midmark 625 has a 650-pound patient weight capacity, which allows you to see more patients in the exam room and to provide those patients with the best care possible.

# Patient-centered care. How does your practice measure up?

More and more of today's patients face a variety of mobility and accessibility challenges, and of course you and your staff want to help patients address those challenges successfully. This can mean re-thinking the workflow so you can focus on bringing care directly to the patient. For example, for many patients the process of being weighed can be uncomfortable, particularly if the scale is in an open hallway, as it often is. With mobility concerns and perhaps waiting for a turn, this pre-exam weighing process can be inefficient, cumbersome and even embarrassing for the patient.

The Midmark 625 Barrier-Free® Examination Table eliminates these challenges with a simple, yet ingenious solution: a digital scale actually built into the body of the exam table.

## Minimize effort. Maximize time. Put patients at ease. The new Midmark 625 with IQscale® makes it all possible.

### A built-in scale for a more effortless exam

The Midmark 625 with IQscale® features a fully integrated and accurate digital scale, which weighs patients quickly, efficiently and discreetly. The handheld controls allow you to run the scale functions (including BMI calculation) with the push of a button and send the data to the EMR. Since patients remain positioned on the exam table, you'll streamline your workflow, save valuable time, and be able to spend more time caring for your patients.

Because weight is recorded electronically and automatically, human error is eliminated; accuracy and privacy are assured. Isn't that a weight off your shoulders?



### All-in-one handheld control

One device does it all. This simple handheld device controls the entire table, including seat height and seatback angle adjustment, scale functions and BMI calculations.



### The conventional weigh-in

For any patient, but especially for those with mobility challenges, the process of being weighed in an ordinary, exposed hallway setting could be an awkward, time-consuming and even embarrassing ordeal.

# Connect patients, practice, records... all in one spot

Consider the future of healthcare: it will mean more connectivity, an increase in integrated data and the expansion of the EMR. In many ways, that future is already here. And the best news is that innovations in equipment are now available to help you and your staff navigate this new technological terrain. The Midmark 625 Barrier-Free® Examination Table offers advanced options like the IQhub®. With the seamlessly integrated IQhub®, gathering data and keeping EMRs up to date will be easier, faster and more accurate. Plus, state-of-the-art connectivity gives the Midmark 625 the ability to better support mandatory Meaningful Use requirements.



## Easy access receptacles and data ports

USB and serial ports on both sides of the table serve as a hub for digital information. Connect directly to vitals devices, so you can take blood pressure, temperature and pulse rate, recording data instantly at the point of care. Your patient's data passes through to your computer to be automatically transferred and stored with the click of a button.



## Hidden connections keep cords neatly out of sight

With the power-source receptacle and data ports located in the rear near floor-level, cords can be kept out of sight for a cleaner look and a safer exam room environment. The rear data ports can be used to connect to laptop computers, creating a data-capturing network directly through the exam table.



## From caregiver to patient. From devices to EMR. Connect your exam room with IQhub®.

### Seamless integration means easier, more efficient examinations

The Midmark 625 with IQscale® and IQhub® streamlines the process of gathering and recording your patients' vitals and keeping EMRs up to date and error free. That means you will be able to focus more of your time and attention on your patients and less time on devices and recordkeeping. Midmark exam tables let you do what you do best – provide excellent care for your patients.

In addition to connecting your vitals, the Midmark 625 with IQhub® lets you connect other digital diagnostic devices at the point of care, including the IQecg® and IQspiro®. The diagnostic data is then ported directly through the IQhub® to your computer with speed and accuracy.





**Powered Multi-Positioning**  
 Custom designed actuators quickly move patients into the proper exam position, allowing you to be more productive. Powered movement of the seat and back gives you the flexibility of positioning patients where you need them.



Barrier-Free® Entry



Supine Position



Lithotomy Position



Low Seated Position




Elevated Seated Position


## New features offer the speed and efficiency you need

### Engineered to maximize proficiency

High performance actuators and new, innovative features like Quick Exam® and Motion Profiling™ give you the freedom to position patients quicker, multi-task when needed, and work more efficiently than any other table on the market.

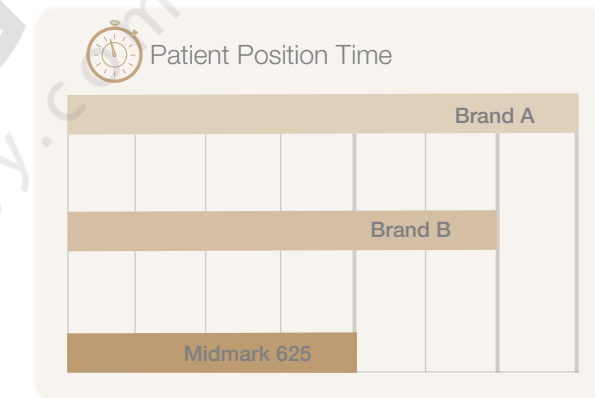


Typical power examination tables take anywhere from 20 to 24 seconds to position patients. The Midmark 625 positions patients in only 14 seconds, saving you valuable time.

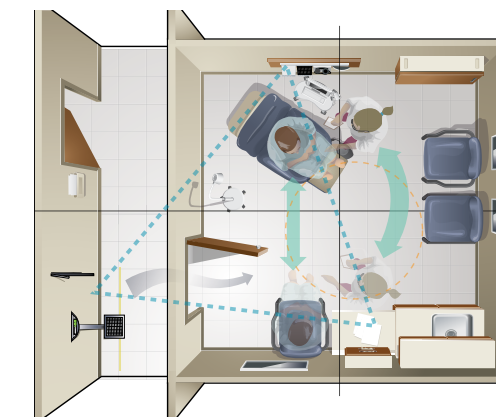

### Traditional workflow: complex and slow

Weighing a patient in a secondary location is time-consuming, cumbersome and perhaps even embarrassing to the patient.

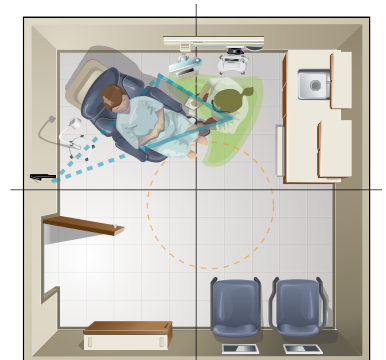

### Today's workflow: simple and streamlined

With its accurate, integrated scale, the Midmark 625 with IQscale® offers a simpler and far more efficient exam protocol. You may also save valuable minutes of time per patient.

## It all comes together right here in the exam room

Meeting new challenges means working more efficiently. And the Midmark 625 Barrier-Free® Examination Table offers everything you need to succeed in the ever-evolving world of healthcare. With full connectivity, generous capacity, a built-in scale and a wide range of convenient features, it all adds up to greater efficiency, accuracy and, ultimately, to a healthier healthcare practice.

# Versatile features and accessories add functionality and efficiency

The Midmark 625 Barrier-Free® Examination Table offers innovative features and accessories designed to give caregivers greater flexibility and help patients feel comfortable and more secure.



## The right amount of storage, organized at the point of care

### All supplies in order and within reach

The patented Exam Assistant® with fold down front gives you unrestricted access to exam supplies, and allows you to get closer to the patient while remaining in an ergonomically correct position. Its unique fold-down front can also be used as a work surface for exam assistance. The larger, overall drawer size allows you to store at least a day's worth of exam supplies. Dividers promote organization for quick supply retrieval.



## Standard Features:



**Hand and Foot Control**  
Standard hand and foot control make patient positioning quick, effective and always within reach.



**Removable Upholstery**  
Allows for color change in seconds or additional access for a thorough cleaning. It's also easy to remove without the use of tools.



**Barrier-Free® Entry**  
Low, flat position reduces the likeliness of patient or caregiver injury.



**Increased Storage**  
Midmark Barrier-Free® tables provide storage for up to four paper rolls, while optional rear storage provides ample room for bulky supplies.



**Premium Comfort System™**  
This patented feature stretches elastomeric material over a seat cutout to enhance patient comfort for a "chair-like" feel on an exam table.



**Electrical Receptacles**  
Hospital grade electrical receptacles on both sides of the table give you easy access to electrical power, freeing up floor space from power cords.



**Motion Profiling™**  
Adjusts speed during backrest movement to mimic the motion of the human body, delivering more natural movement and reducing the time needed to position the patient.



**Quick Exam®**  
The Quick Exam® feature allows you to elevate the table to the exam position with one touch of the control.

## Available Accessories:

### Wireless Controls

Remove cords and eliminate obstacles. Enjoy the freedom of locating the hand and foot control anywhere in the room. Low battery audible indicator provides early warning of approaching battery replacement.



Wireless Hand Control



Wireless Foot Control




**IQvitals® Table Mount**  
Convenient table mount for the Midmark IQvitals® vital signs device, designed to fit all Midmark and Ritter exam tables.



### Adjustable Arm System

Offers patients stable support for accessing the exam table and a sense of security during transfer. Arms rotate and lock into one of three positions and are easily removable.

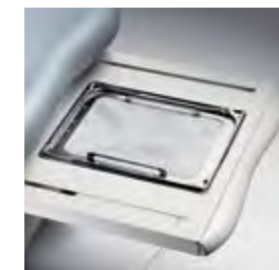

**Urology Drain Pan**  
Ideal for urology procedures or other situations in which removal of fluid and debris is necessary.



### Assist Handles

Provide patients a secure means of access to the exam table or to reposition themselves during an exam. Easily removed without tools, assist handles offer greater patient access.



**Supply Assistant**  
Provides extra storage for bulky exam supplies. Keeps drapes, gowns and similar items within easy reach.



### Articulating Knee Crutch

A more ergonomic and comfortable option compared to standard stirrups. The ball joint design offers ultimate flexibility in patient positioning. Available in a non-articulating option.



**Ritter 253 LED Exam Light**  
Easy access controls, high intensity light and an adjustable focal spot makes the 253 the perfect complement to any procedure room. Available in table mount and mobile options.

## Specifications

Patient Weight Capacity:  
650 lb (295 kg)

Height:  
Minimum 18" (46 cm)  
Maximum 37" (94 cm)

Premium Upholstery Dimensions (top):  
28" W x 59.5" L  
(71 cm x 151 cm)

Length with footrest extended:  
75.5" (190 cm)

Stirrups are standard

Pelvic Tilt: 7°


Paper Roll Holder holds 21" x 3.5"  
dia. paper roll (53 cm x 9 cm)

Electrical Receptacle:  
(115 VAC, 5 A Max.)

Hospital grade, electrical  
receptacles on left and right sides

Table Weight:  
433 lb (197 kg)

Shipping Weight:  
505 lb (229 kg)

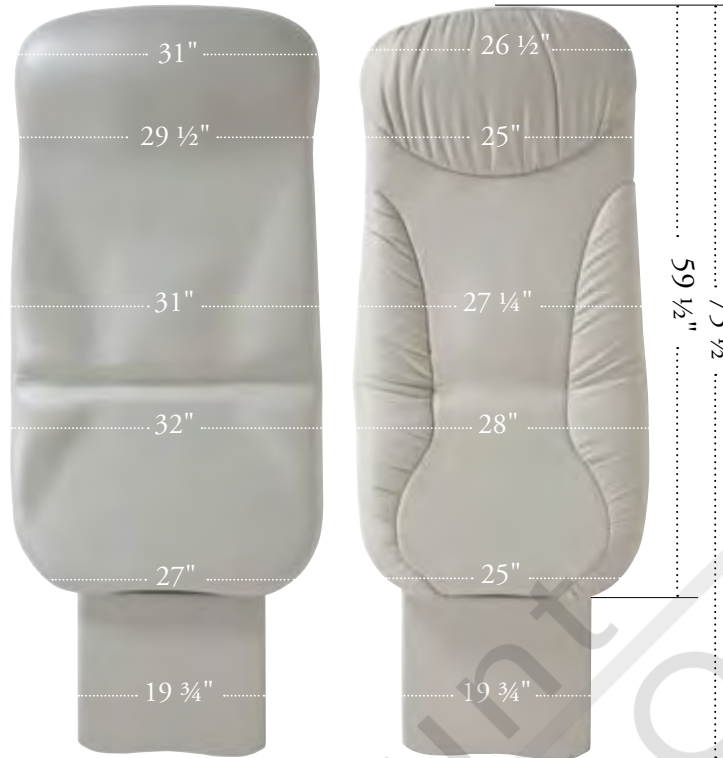

## Electrical Ratings

625-001 – 115 VAC, 50/60 Hz, 12 A

625-003 – 115 VAC, 50/60 Hz, 12 A

625-004 – 115 VAC, 50/60 Hz, 12 A

## Safety Standards

UL 60601-1, CAN/CSA-C22.2

No.601.1-M90, EN 60601-1

NFPA 99

## Options

**Ultra-Premium Top:**  
28" W x 59.5" L (71 cm x 151 cm)

**Wide Premium Top:**  
32" W x 59.5" L (81 cm x 151 cm)  
Paper Roll Holder holds  
24" x 3.5" dia. paper roll  
(61 cm x 9 cm)

**IQscale®:**  
Scale Capacity: 650 lb (294.8 kg)  
Scale Resolution: 0.2 lb (0.1 kg)

Scale Accuracy:  
30 - 100 lb  
(13.6 - 45.4 kg) ± 0.2 lb (0.1 kg)

100 - 300 lb  
(45.4 - 136.1 kg) ± 0.4 lb (0.2 kg)

300 - 500 lb  
(136.1 - 226.8 kg) ± 0.6 lb (0.3 kg)

500 - 650 lb  
(226.8 - 294.8 kg) ± 0.8 lb (0.4 kg)

**IQhub®:**  
USB ports on left and right sides,  
and rear of table

USB speeds up to 480 Mbit/sec

Vitals device ports on left and right  
sides, plus rear of table

## Accessories

**Supply Assistant:**  
22.4" W x 17.9" L x 12.1" H  
(57 cm x 46 cm x 31 cm)

**Adjustable Arm System:**  
Right and left side arms included

**Assist Handles:**  
Right and left side assist  
handles included

**Knee Crutches:**  
articulating and standard

## PC Work Shelf

**Flat Upholstery:**  
Stitched, without pillow

## Urology Pan With Hose

## Stainless Steel Pan

## Premium Solid Upholstery Colors

Blueberry-235

Cashmere Blue-239

Silver Sage-236

Terra Cotta-240

Fossil Grey-238

Tea Green-237

Black-312

## Premium Pattern Upholstery Colors

Belagio-247

Capri-251

Siena-248

Treviso-252

Firenze-250

Tuscany-249

## Ultra-Premium Solid Upholstery Colors

Ultra-Pacific Blue-241

Ultra-Sky-245

Ultra-Spruce-242

Ultra-Canyon-246

Ultra-Arctic-244

Ultra-Olivine-243

Ultra-Raven Wing-228

CARB 93120.2 Phase 2 Compliant

Midmark is an ISO 13485 and ISO 9001 Certified Company  
For more information or a demonstration, contact your Midmark  
dealer or call: 1-800-MIDMARK Fax: 1-800-365-8631  
Outside the U.S.A. call: 1-937-526-3662 Fax: 1-937-526-8392  
or visit our website at midmark.com

© 2009 Midmark Corporation  
60 Vista Drive P.O. Box 286 Versailles, Ohio 45380-0286  
Products subject to improvement changes without notice  
Litho in U.S.A. 007-0828-00 Rev. D1 (12/13)



Because we care.