




Barrier-Free® Power Examination Table

623



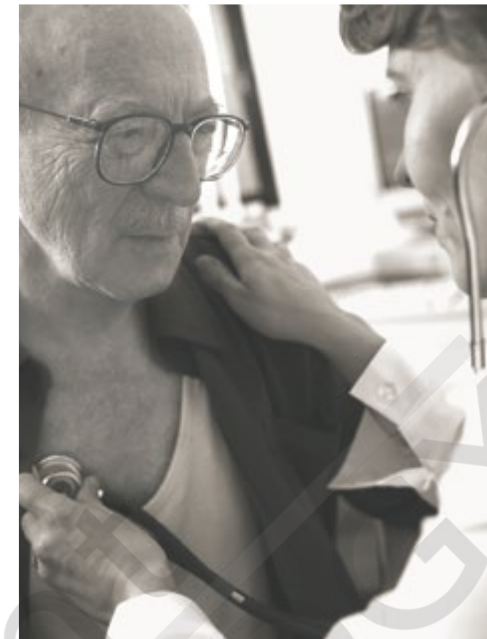

Break down the barriers®

The face of healthcare is changing - and so are your patients.

The percentage of patients you see who are elderly, disabled, obese or have other mobility concerns is steadily on the rise. To keep up with this shifting demographic, you must make changes in healthcare delivery. You can start with a simple step: equipping your exam rooms with products that accommodate all patients.

Your patients come in all shapes and sizes.

We typically think of the elderly or patients with disabilities that are obvious as needing help getting onto or off of the exam table – but that's not always the case. Take, for instance, children. With a low height of 18 inches from the floor, Midmark Barrier-Free® exam tables make accessing the exam table a snap regardless of height or physical condition.



Accessibility for all patients, by design

The patients you see are all different. But the one thing they may have in common is difficulty accessing the exam table. Easy access to the exam table can decrease strain on patients and staff, and ensure that you can position all patients effectively for a more thorough and accurate exam. Break down the barriers in the exam room with a Midmark Barrier-Free® exam table, and provide all of your patients with equal care regardless of physical condition.



Whether transferring a wheelchair patient or assisting an expectant mom, one thing is clear: the Midmark 623 means that positioning patients onto the exam table will be easier, more efficient, and much safer for you and your staff. With its 18- to 37-inch adjustable table height, lifting is less stressful, more controlled, and the patient's dignity is preserved.

Comfort, safety and efficient patient care

To meet the needs of all patients, including those who are disabled or those with mobility concerns, as well as the physicians who care for them, Midmark Barrier-Free® examination tables are designed to decrease the chance of injury and reduce patient anxiety.

"Decreasing the need for assistance gives the patient an improved sense of confidence."

By listening to the needs of physicians and their patients, Midmark has developed products that deliver higher levels of exam room efficiency and comfort to the patient, the physician and the staff. Barrier-Free® exam tables also reduce the likelihood of distress and injury to patients who may have difficulty accessing a fixed height exam table, such as the elderly, obese or even expectant mothers.



Patient positioning made easy

Custom designed actuators quickly move patients into the proper exam position, increasing efficiency. Powered movement of the seat and back gives you the flexibility of positioning patients where you need them.



Standard Features:



Multi-function Foot & Hand Control

The foot control moves the seat and back section independently or simultaneously for the most efficient patient movement. The 623 hand control gives an alternative method of positioning patients.



Accessible Controls

Upholstery and drawer heater controls are easily accessible for you and your staff.



Removable Upholstery

Allows for color change in seconds or additional access for a thorough cleaning. It's also easy to remove without the need for tools.



Premium Comfort System™

This patented feature stretches elastomeric material over a seat cutout to enhance patient comfort for a "chair-like" feel on an exam table.

Available Options:



Ultra-Premium Upholstery
Offer sophistication and comfort like never before with a supple, soft texture resembling the finest calfskin leather. And, as an added bonus, it's PVC-free. (See the back page for color options.)



Heated Upholstery
Heated seat and lumbar improve patient comfort and reduce patient anxiety.

Accessories when you want to add versatility and expand functionality

The Midmark 623 Barrier-Free® Examination Table offers innovative accessories designed to give caregivers greater flexibility and help patients feel comfortable and more secure.



Urology Drain Pan
Ideal for urology procedures or other situations in which removal of fluid and debris is necessary.




Stainless Steel Treatment Pan
Provides you with a safe place to store instruments and supplies at the point of care.



Adjustable Arm System
Offers patients stable support for accessing the exam table and a sense of security during transfer. Arms rotate and lock into one of three positions and are easily removable.




Increased Storage
Midmark Barrier-Free® tables provide storage for up to four paper rolls, while optional rear storage provides ample room for bulky supplies.

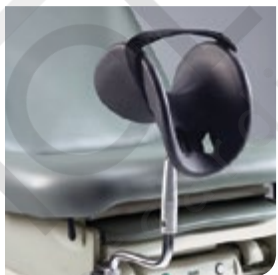

Standard Knee Crutches
Swivel to provide leg support for extended procedures where fatigue is a factor or when the patient lacks sufficient leg strength or control.



Supply Assistant
Provides extra storage for bulky exam supplies. Keeps drapes, gowns and similar items within easy reach.

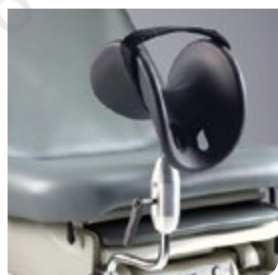

Articulating Knee Crutch
A more ergonomic and comfortable option compared to standard stirrups. The ball joint design offers ultimate flexibility in patient positioning. Available in a non-articulating option.



Assist Handles
Provide patients a secure means of access to the exam table or to reposition themselves during an exam. Easily removed without tools, assist handles offer greater patient access.



PC Work Shelf
Provides an adjustable height workspace for a tablet PC, portfolio, or for exam instruments and supplies.



IQvitals® Table Mount
Convenient IQvitals® swing arm table mount positions vitals at the side of the exam table to support efficient workflow and keeps cords secure and out of the way.

Specifications

Dimensions (top):
 28.5" W x 59.5" L
 (72.4 x 151.1 cm)
 Length with footrest extended:
 75.5" (191.8 cm)
 Height: Minimum 18" (45.7 cm)
 Maximum 37" (94.0 cm)

Patient Weight Capacity:
 400 lb (181.4 kg)

Paper Roll Holder:
 Holds 21" x 3.5" paper roll
 (53.3 x 8.9 cm)

Electrical Receptacle:
 (115 VAC 5 Amps Max.) Duplex,
 hospital grade, electrical receptacle
 on left side

Rear storage pod:
 (includes removable divider)
 17.8" W x 7.0" L x 9.8" H
 (45.2 x 17.8 x 24.8 cm)

Stirrups are standard.

Electrical Requirements – 623

115 VAC, 50/60 Hz, 13.5 Amps
 (with heaters – 008)

Safety Standards

UL 60601-1 and
 CAN/CSA-C22.2 No.601.1-M90

Options

Premium/Ultra-Premium Top
 28.5" W x 59.5" L
 (72.4 x 151.1 cm)

Pelvic Tilt (5 degrees)

Heated Upholstery Top
 Temperature Range:
 98°F (37°C) to 105°F (41°C)

Heated Drawer (35 watt)
 Temperature: 98°F (37°C)

Accessories

Adjustable Arm System
 (right and left side arms included)
 Factory installed only

Articulating Knee Crutches

Standard Knee Crutches

PC Work Surface

Supply Assistant

LED Exam Light

Urology Pan

Stainless Steel Pan

Ritter 253 LED Exam Light



Ritter 253 LED Exam Light

Easy access controls, high intensity light and an adjustable focal spot makes the 253 the perfect complement to any exam room. Available in table mount and mobile options.

Premium Solid Upholstery Colors

Blueberry-235

Cashmere Blue-239

Silver Sage-236

Terra Cotta-240

Fossil Grey-238

Tea Green-237

Black-312

Premium Pattern Upholstery Colors

Belagio-247

Capri-251

Siena-248

Treviso-252

Firenze-250

Tuscany-249

Ultra-Premium Solid Upholstery Colors

Ultra-Pacific Blue-241

Ultra-Sky-245

Ultra-Spruce-242

Ultra-Canyon-246

Ultra-Arctic-244

Ultra-Olivine-243

Ultra-Raven Wing-228

CARB 93120.2 Phase 2 Compliant

Midmark is an ISO 13485 and ISO 9001 Certified Company.
 For more information or a demonstration, contact your Midmark dealer or call: 1-800-MIDMARK Fax: 1-800-365-8631
 Outside the U.S.A. call: 1-937-526-3662 Fax: 1-937-526-8392
 or visit our website at midmark.com

© 2005 Midmark Corporation
 1700 S. Patterson Blvd., Suite 400, Dayton, Ohio 45409
 Products subject to improvement changes without notice
 Litho in U.S.A. 007-0590-00 Rev. F1 (3/14)

