

# T1

## COMPACT TRANSPORT MONITOR

### Key Features


- 5-inch brilliant touchscreen display and robust monitoring capability in compact 2.0 lb (0.9 kg) package
- Up to 4 waveforms
- Standard features include 3 or 5-lead ECG, SpO<sub>2</sub>, NIBP, respiration, temperature, dual invasive pressure
- Choice of Masimo SET® SpO<sub>2</sub> or Nellcor® OxiMax® SpO<sub>2</sub>
- Optional features include
  - ST and arrhythmia analysis
  - Microstream® or DPM CO<sub>2</sub>
- Connections for external devices
- Extensive data storage capability for trend data, alarms, events, NIBP measurements and up to 48 hours of full disclosure
- Docking Station connection for battery charging and printing real-time and historical data
- Ergonomic carrying handle and multiple mounting options offer unprecedented mobility
- Transport Module Rack offers expanded functionality



## your compact mobile assistant

High quality medical care deserves a transport solution that not only is easy to use and easy to carry, but also provides the advanced monitoring capabilities needed within your hospital environment. With its brilliant touchscreen display and carefully designed features, the Mindray T1 has been designed to meet the unique challenges of monitoring patients during transport.

# T1 COMPACT TRANSPORT MONITOR

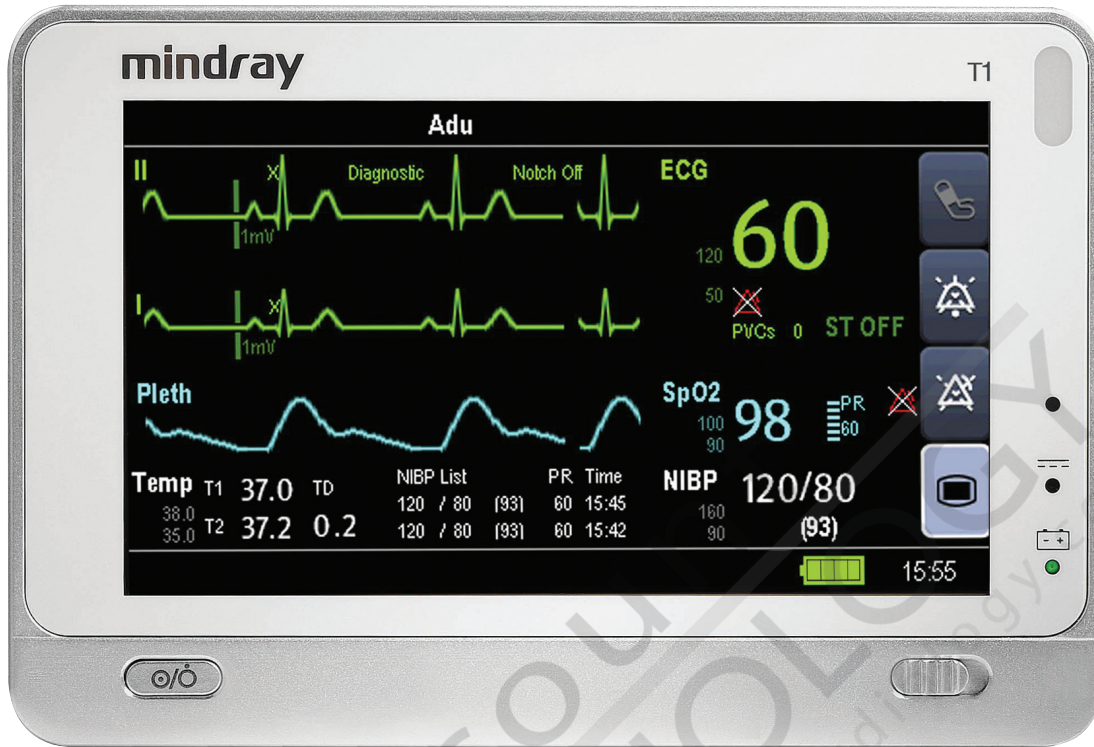

5-inch high-resolution touchscreen display and advanced monitoring in a compact 2.0 lb (0.9 kg) package.

*actual size*



Add functionality with the Transport Module Rack



For flexibility, the T1 can be easily mounted to a bed rail or pole



Docking Station offers connection to a networked laser printer and charges the battery

Mindray DS USA, Inc.  
800 MacArthur Blvd., Mahwah, NJ 07430  
Tel: 800.288.2121 Tel: 201.995.8000 Fax: 800.926.4275 [www.mindray.com](http://www.mindray.com)

©2014 Mindray DS USA, Inc.  
Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 5/14 P/N: 0002-08-30023 Rev F

**mindray**  
NORTH AMERICA