

Accutorr® V

VITAL SIGNS MONITOR

Features and Benefits

- Large bright LEDs are easily read from any vantage point, under various lighting conditions – even at a distance
- Dedicated buttons for all primary functions provide incomparable ease of use
- 3.4" high resolution color LCD displays trend data, plethysmograph, alarms, events and menus
- Every unit includes a lithium ion battery, a built in recorder and a pulse strength indicator
- Many versatile options include a bar code scanner, Masimo SET® or Nellcor® OxiMax® SpO₂; Mindray SmarTemp™ or Welch Allyn® SureTemp® Plus temperature
- Patient size, inflation pressure, interval time, alarms, audio level, date and temperature format are all easily configured
- Takes blood pressure just like you do, with a gentle linear bleed; readings are easily confirmed with a stethoscope
- Unique, flexible database makes up to 1,200 readings instantly accessible by patient ID, if desired
- Supports connectivity to e-Gateway for data acquisition and storage to your electronic medical record
- Supports nurse call and industry-standard RS-232 communications for connection to other medical devices via our proprietary DIAP protocol



confident | capable | comfort

The Accutorr V vital signs monitor boasts a sleek design with more features and more capability. It offers fast, direct operation with the comfort and ease that you would expect from an Accutorr. It's all there. Every unit comes equipped with a high resolution color LCD, recorder, lithium ion battery and all the accessories you need to begin monitoring. Options include your choice of SpO₂ technology, including Masimo SET® or Nellcor® OxiMax® SpO₂ and temperature. Plus, by using e-Gateway or the Accutorr CS charting solution, you can conveniently connect to your facility's EMR.


Accutorr V offers proven Accutorr performance with our most advanced NIBP technology, designed specifically for ease of use – all in a lightweight, portable package. If maximum flexibility, high value, and reliability match your vital signs needs, then look no further than Accutorr V.

Accutorr® V

The intuitive software of the Accutorr CS eliminates manual charting and allows for the direct import of vital signs data into a patient's Electronic Medical Record (EMR). It links the point-of-care directly to the nurse's terminal thus improving workflow and reporting accuracy. Now you can create a new standard of care with the Accutorr V vital signs monitor and Accutorr CS charting solution.



Automatic trend graphs and tables can be viewed on screen or sent to your laser printer.

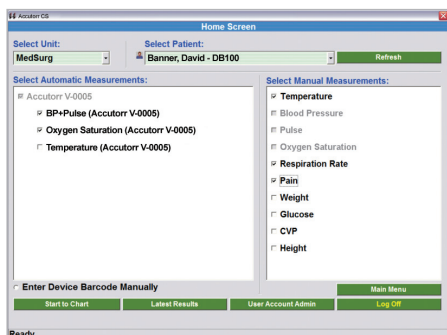

REV ID: 20011101170906								Adult
Date	Time	Sys	Dia	Map	SpO2	PR	Temp	
01-01	23:06	123	/ 80	/ 94	98	73	98	
01-01	23:05	114	/ 76	/ 88	97	76	98	
01-01	23:04	121	/ 75	/ 87	98	79	98	
01-01	23:03	137	/ 74	/ 95	97	77	97	
01-01	23:02	125	/ 67	/ 87	98	79	98	
							01-01-2011	17:11:38

Color LCD supports trend or waveform display

SmarTemp™



Bar code ready – plug in a bar code scanner and track vitals by patient ID.



Welch Allyn® SureTemp® Plus

Temperature – at the speed you want and the accuracy you demand – whether using the optional SmarTemp™ or Welch Allyn® SureTemp® Plus thermometer.

Mindray DS USA, Inc.
 800 MacArthur Blvd., Mahwah, NJ 07430
 Tel: 1.800.288.2121 Tel: 201.995.8000 Fax: 1.800.926.4275 www.na.mindray.com

Mindray™ and SmarTemp™ are trademarks of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. Accutorr® is a trademark of Mindray DS USA, Inc. and is registered in the U.S. and other countries. Masimo SET® is a U.S. registered trademark of Masimo Corp. Nellcor® and Oximax® are registered trademarks of Nellcor Puritan Bennett, Inc. Welch Allyn® and SureTemp® are registered trademarks of Welch Allyn, Inc. Datascope is a trademark or registered trademark of Datascope Corp. or its subsidiaries and is used under license by Mindray DS USA, Inc. © 2012 Mindray DS USA, Inc. Subject to change. P/N: 0002-08-6353 Rev J

