


Summary of Features and Benefits

- 12.1" TFT display, configurable up to 8 waveforms, enables maximum data analysis of your most critical parameters
- Features include 3/5-lead ECG, non-invasive blood pressure, Nellcor® OxiMax® SpO₂, pulse rate, respiration, dual-channel invasive blood pressure, dual-channel temperature, lithium ion battery, and integrated recorder
- Optional features include DPM Sidestream CO₂, and arrhythmia and ST analysis
- 96-hours of graphic and list trends, NIBP list including 800 measurements, and alarm event recall for 70 alarm episodes for comprehensive historical review of all parameters
- 40-seconds of full disclosure waveform recall for detailed review of physiological events
- Compact flash back-up, included on every device, protects your patient's data should a sudden loss of power occur
- Built-in, 3-trace recorder for printing real-time and historical data

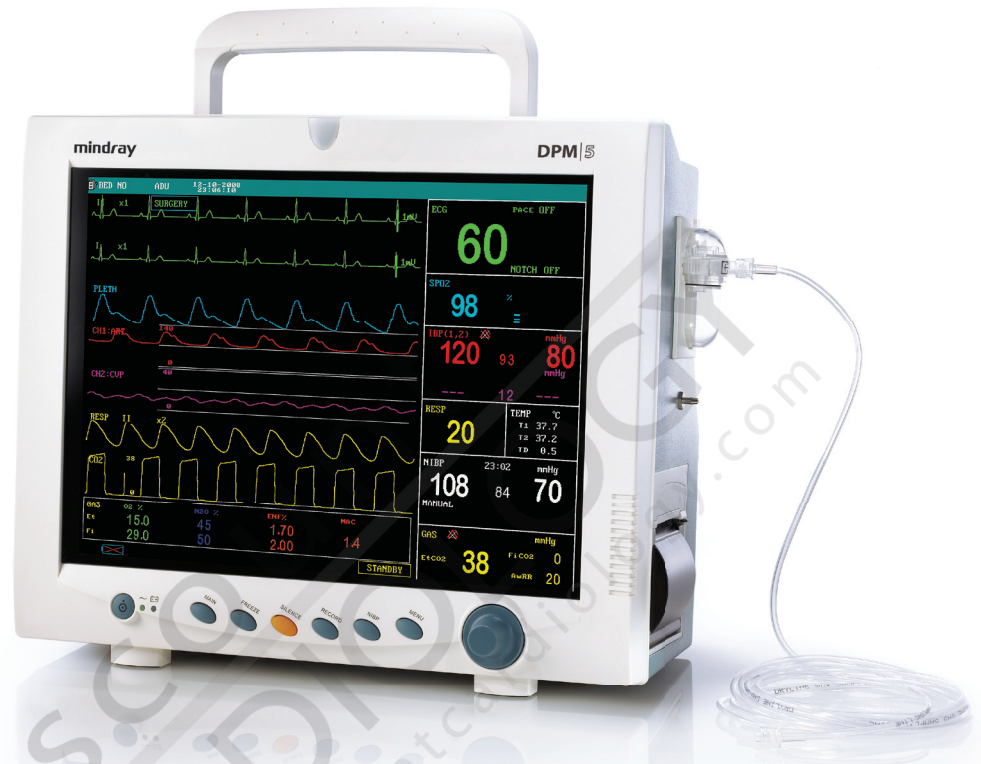


Highly capable, yet conveniently flexible.


DPM5 is a robust patient monitor, small enough to move around, yet flexible enough to adapt to your particular needs. With optional features such as ETCO₂ measurement and arrhythmia and ST analysis, DPM5 is configurable to your liking. Quick-action keys on the front panel make this device especially easy to learn, which means a shorter learning curve for new or changing staff. Rich in function, conveniently adaptable and pleasingly affordable, DPM5 may be just the monitor for your department.

DPM™ | 5

PATIENT MONITOR

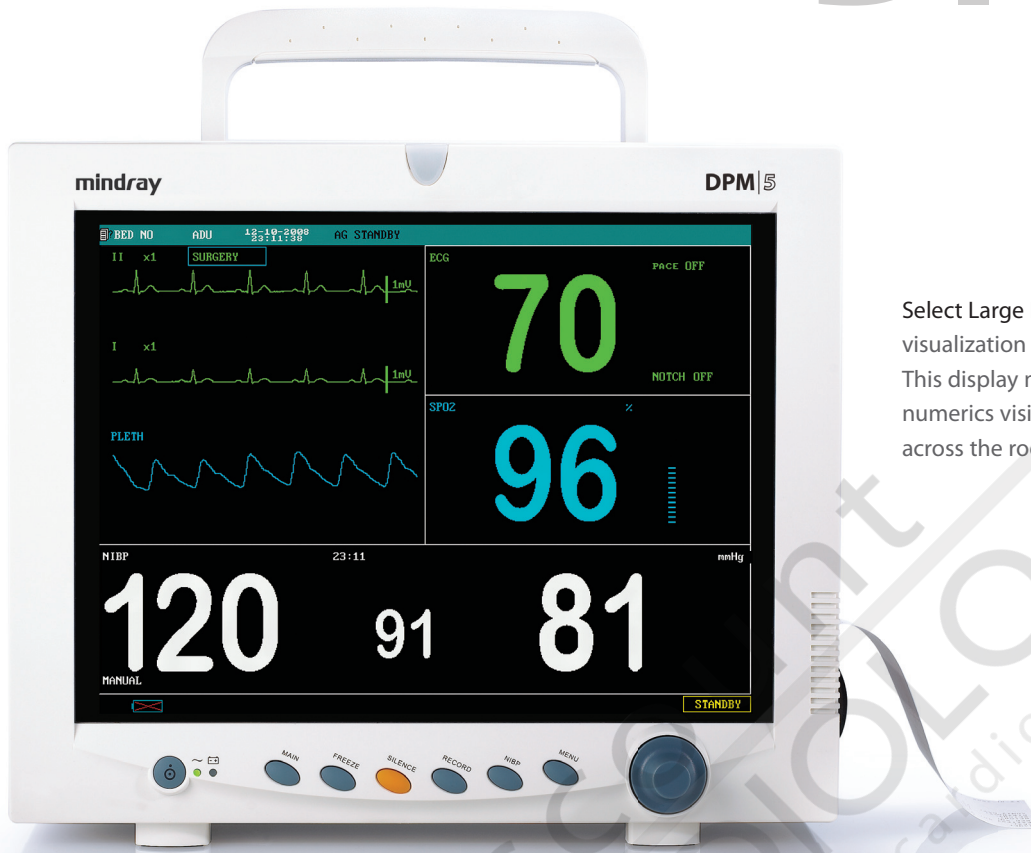

Select Large Numerics Mode for optimal visualization of your most important data. This display mode features bright, sizable numerics visible from across the bed or across the room.

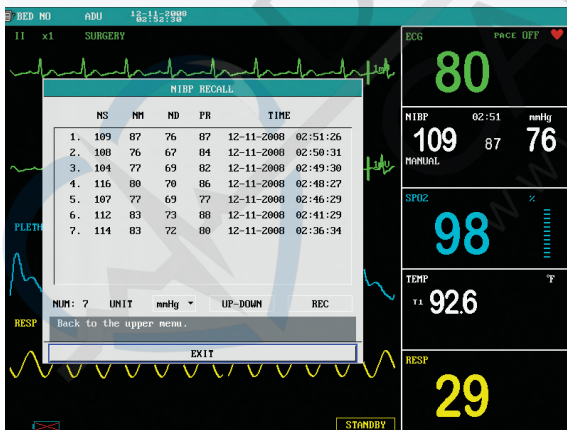

Review your patient's most recent blood pressure measurements using the NIBP recall function, accessible using the Menu quick-action key. At the touch of a button, obtain a dedicated printout of your patient's NIBP measurements.



Choose from a series of mounting schemes to create your ideal monitoring environment.

Mindray DS USA, Inc.
 800 MacArthur Blvd., Mahwah, NJ 07430
 Tel: 1 800 288 2121 Tel: 201 995 8000 Fax: 1 800 926 4275

www.na.mindray.com

mindray™
 NORTH AMERICA