


CMS Viewer

Features and Benefits

- A software-based solution providing clinicians with remote access to critical patient data on the DPM Central Station network
- Review near real-time streaming of patient data from virtually anywhere through the hospital VPN
- 16 simultaneous users
- No dedicated server required
- Physiological alarms on target monitor are clearly identified on the remote display based upon alarm priorities
- Through a secure hospital provided VPN connection, the CMS Viewer provides direct access to all patient waveforms and parameters captured on the DPM Central Station Network
- An additional option is the Watchdog Assembly. This feature is a single use license, sold separately, which expands the capability of the CMS Viewer to review trends, events and alarms for current and discharged patients.

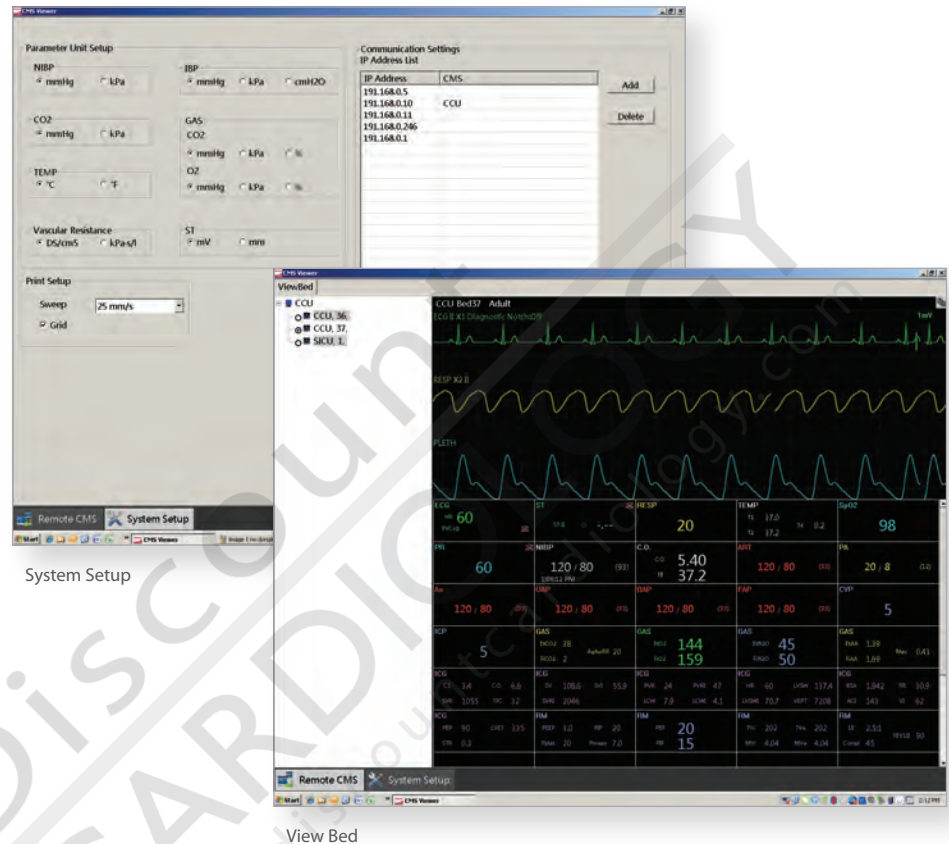

Remote access web viewer for the DPM™ Central Station

The Mindray CMS Viewer expands the range and scope of the DPM Central Station, allowing critical access to timely information within the hospital network or through a VPN connection across the world wide web. Authorized users install the CMS Viewer software, access the hospital network and have access to all authorized patients on the DPM Central Station network.

CMS Viewer with Watchdog Assembly

Watchdog Assembly

This single user license adds the following capabilities:

- Review of waveforms, trends, alarms, and events of current and discharged DPM CS patients
- Print and store waveform data

To include:

- View Bed – Near real-time viewing of all patients
- Waveform Review – 240 hours in 20 minute compressed format viewing up to 12 leads
- Trend Review – graphic or tabular
- Event Review – By event or priority type
- ST Review – View and print historical ST Segment
- C.O. Review – Review as many as 720 past measurement curves
- Patient Management – Patient Demographics


System Requirements

CMS Viewer PC Platform Minimum

- Windows® 2000, Windows® XP
- Intel Pentium IV 2.0 GHz
- 512 MB of RAM
- 40GB hard disk, 7200 rpm SATA
- Ethernet 802.3, RJ45 port

Display Requirements Minimum

- 17" LCD minimum 1280 x 1024 (dpi)
- Length to width ratio of 4:3 or 5:4 with CE marking

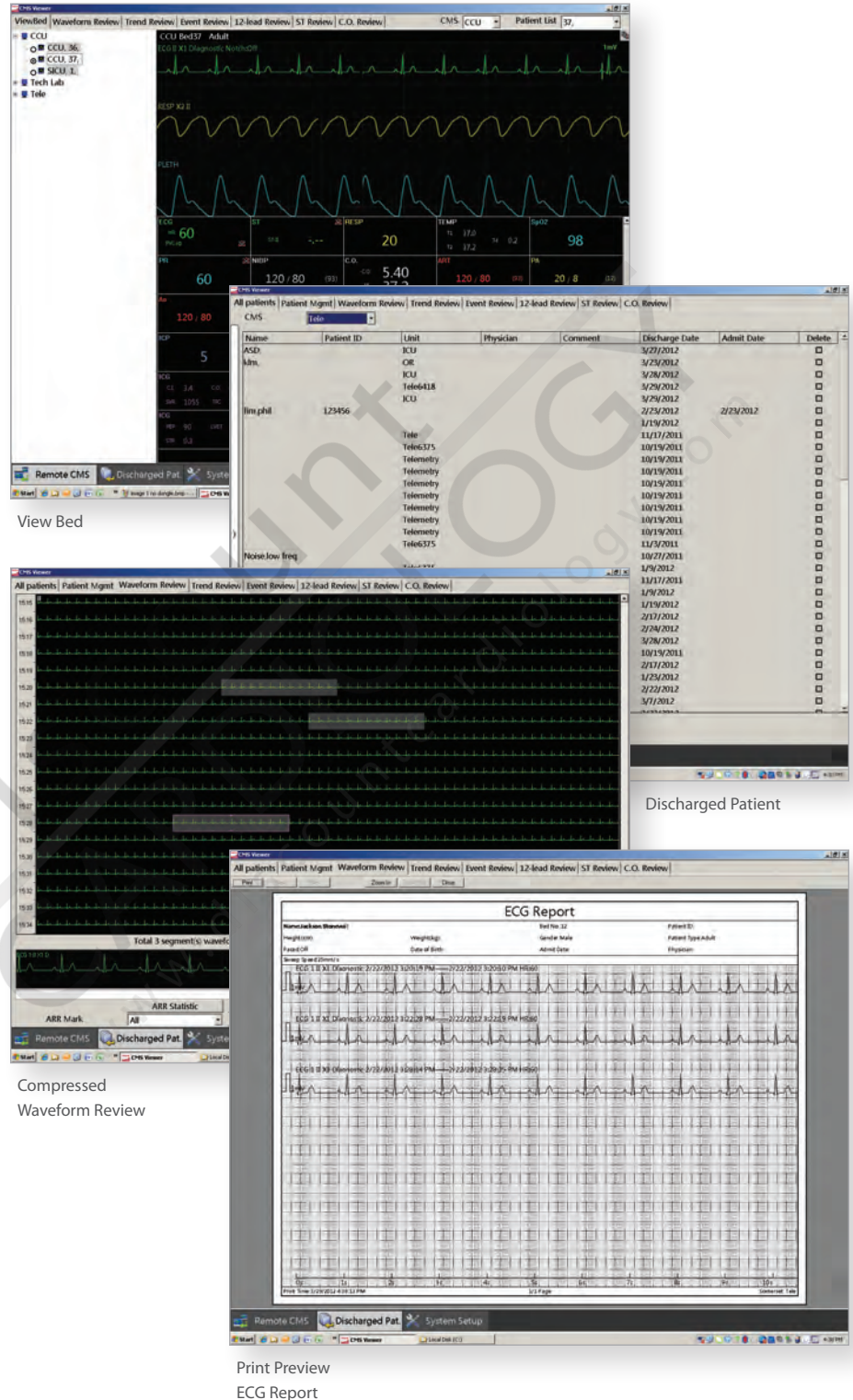

Mindray DS USA, Inc.
 800 MacArthur Blvd., Mahwah, NJ 07430
 Tel: 1.800.288.2121 Tel: 201.995.8000 Fax: 1.800.926.4275 www.na.mindray.com

12-Lead Review Analysis tab not available in U.S.
 DPM™ is a trademark of Mindray DS USA, Inc. Mindray™ is a trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd.
 Windows® is a registered trademark of Microsoft Corporation in the United States and other countries.
 ©2012 Mindray DS USA, Inc. Subject to change. P/N: 0002-08-37010 Rev C