

## Summary of Features and Benefits

- Weighing less than 8 lbs., its lightweight design provides a highly mobile, vital signs solution
- Smart, stable design anchors DPM3 firmly in place when activating buttons, keeping your second hand free to perform other tasks
- Standard parameters include non-invasive blood pressure, SpO<sub>2</sub>, and pulse rate
- Optional parameters include Masimo SET® or DPM SpO<sub>2</sub> and temperature
- Sizable, bright LEDs provide optimal visibility, whether across the bed or across the room
- 3.4" high resolution LCD for display of plethysmograph, trended parameter data, alarm events and menus
- With SmarTemp™ integrated temperature perform predictive oral, axillary and rectal measurements in less than 12 seconds from calibration
- Data storage is easy and abundant, storing up to 100 unique patient IDs and 1,200 sets of data for each parameter. Plus, all data is transferable to a PC for complete review and printing.
- Built-in recorder comes standard for printing of real-time and historical data
- Automatically admit patients using the optional bar code scanner. The process is quick and efficient, and can aid in minimizing risk of common transcription errors.



The unusually affordable and refreshingly intuitive vital signs monitor.

DPM3 is a compact, easy-to-use vital signs monitor designed to satisfy your most basic monitoring needs. Suitable across adult, pediatric and neonatal patient populations, DPM3 is well-adapted for any department where quick, reliable physiological measurements are required. And, since it's unusually affordable, you can consider DPM3 for any area of the facility in which vital signs may be needed – whether on a regular or sporadic basis. Equipped with lithium ion battery technology that provides over 10 hours of continuous run-time, this powerful little monitor can keep up with you on even your busiest shifts. With built-in carrying handle, rolling stand, wall, IV pole and bed mounting options, accompanied by a robust, lightweight design, DPM3 is ready to serve as your everyday vital signs sidekick.

# DPM™ | 3

VITAL SIGNS MONITOR

Alarm light indicates three levels of user-adjustable alarms

Vivid, 0.8" LEDs prominently display NIBP measurement

Built-in thermal recorder comes standard on every device


SmarTemp™ thermometer with automatic activation when probe is removed

Patient size indicator confirms whether in adult, pediatric or neonate mode

Temperature indicator confirms whether in predictive or monitor mode, as well as selected measurement site

One-touch buttons for access to all major functions

High resolution, 3.4" LCD displays plethysmograph, trended data and alarms

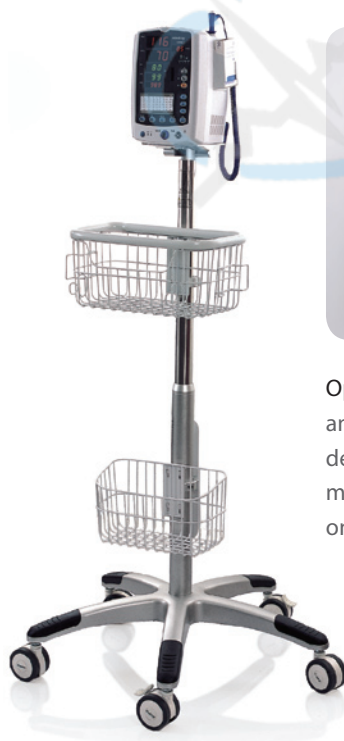

Choose from 2 temperature measurement modes, predictive or monitor, depending on your clinical needs.

Optional rolling stand, with ample storage space and sturdy design, provides a convenient means of performing vital signs on the go.

Mindray DS USA, Inc.  
800 MacArthur Blvd., Mahwah, NJ 07430  
Tel: 1 800 288 2121 Tel: 201 995 8000 Fax: 1 800 926 4275

[www.mindray.com](http://www.mindray.com)

**mindray™**  
NORTH AMERICA