



Critical Values

Speed the delivery of urgent care

Get critical information at the point of care:

- Quickly identify patients that need urgent care
- Shorten the time from discovery to treatment
- Supports The Joint Commission guidance

Philips PageWriter TC70 Cardiograph makes it simpler than ever to align with The Joint Commission's National Patient Safety Goals to improve the timeliness of reporting critical test results and their receipt by a responsible licensed caregiver.

When appropriate, the PageWriter TC70 automatically displays Critical Value summary statements boldly on screen and on printed ECG reports to advise caregivers of the need for urgent care. Clinicians can quickly take action and reduce the time from the discovery of a critical cardiac event to intervention, such as balloon angioplasty.

The PageWriter TC70 summarizes 30 interpretive statements into four, easily understood, Critical Values:

- Acute MI
- Acute Ischemia
- Complete Heart Block
- Very High Heart Rate


This user-friendly system speeds the care process from downloading orders, taking an ECG, and identifying critical test results, to enable clinicians to quickly provide needed treatment.

PHILIPS

Identify critical test results quickly and easily

Acute MI

Anterior infarct, acute
 Probable anterior infarct, acute
 Anterior infarct, possibly acute
 Extensive anterior infarct, acute
 Extensive anterior infarct, possibly acute
 Inferior infarct, acute
 Probable inferior infarct, acute
 Inferior infarct, possibly acute
 Inferoposterior infarct, acute
 Lateral infarct, acute
 Probable lateral infarct, acute
 Lateral infarct, possibly acute
 Inferolateral infarct, acute
 Anterolateral infarct, acute
 Probable anterolateral infarct, acute
 Anterolateral infarct, possibly acute
 Anteroseptal infarct, acute
 Probable anteroseptal infarct, acute
 Anteroseptal infarct, possibly acute
 Posterior infarct, acute
 Probable posterior infarct, acute
 Right ventricular infarct, acute
 Probable right ventricular infarct, acute

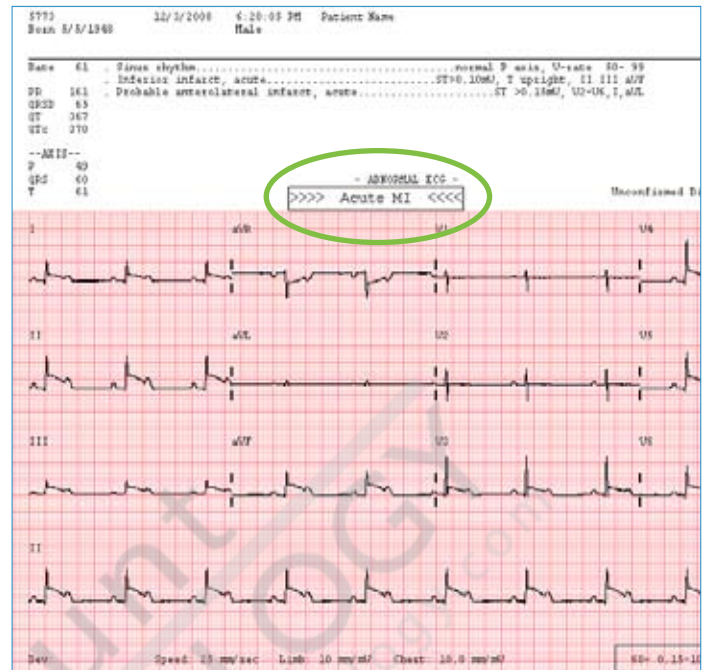

Critical Value statements are boldly displayed on screen and in printed ECG reports

Acute Ischemia

Repolarization abnormality, severe global ischemia

Complete Heart Block

AV block, complete (third-degree)
 Complete AV block with wide QRS complex
 Atrial-flutter/fibrillation with complete AV block

Very High Heart Rate

Extreme tachycardia
 Wide-QRS tachycardia
 Extreme tachycardia with wide complex



© 2008 Koninklijke Philips Electronics N.V.
 All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Philips Healthcare is part of Royal Philips Electronics

www.philips.com/healthcare
healthcare@philips.com
 fax: +31 40 27 64 887

Printed in The Netherlands
 4522 962 39051/860 *SEP 2008

Philips Healthcare
 Global Information Center
 P.O. Box 1286
 5602 BG Eindhoven
 The Netherlands