


Capnocheck® II Hand-Held Capnograph/Oximeter



The BCI® Capnocheck® II Hand-Held Capnograph/Oximeter is a full-featured, highly portable unit available in your choice of configurations, delivering capnography only or capnography combined with oximetry. Ideally suited for use in and out of the hospital, its back-lit display assures easy viewing of wave-forms, numeric values and on-screen trending data. The battery-operated design with an optional battery bypass offers the flexibility to accommodate both stationary and transport applications. Plus, an infrared link to an external printer provides quick, trouble-free downloading.

For more information visit smiths-medical.com/bci.

PATIENT MONITORING

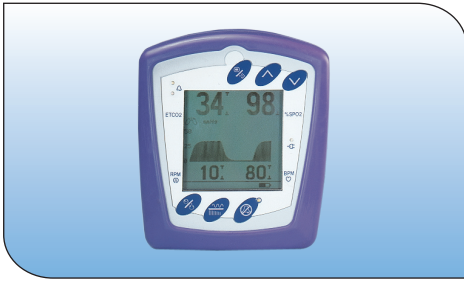

Capnocheck® II Hand-Held Capnograph/Oximeter

- Measures ETCO₂, inspired CO₂, respiration rate, SpO₂ and heart rate.
- Provides waveforms, numeric values and on-screen trending.
- Easy-to-read EL back-lit LCD display.
- Adjustable alarm limits, as well as visual and auditory alarms.
- Sidestream technology accommodates intubated and non-intubated patients.
- Lightweight for easy transport.
- Infrared link to an external printer assures trouble-free download.
- Monitor patients ranging from pediatric to adult.
- Genuine BCI® reusable and disposable sensors available.

Technical Specifications

DISPLAYS, INDICATORS AND KEYS

- **Display**
LCD, with EL back light
- **Indicators**
High and medium/low priority alarm, Alarms silenced, and AC power
- **Displayed Trends**
Graphical display for ETCO₂, respiration rate, % SpO₂, or pulse rate; user selectable scale
- **Displayed Waveforms**
CO₂, waveform or oximetry plethysmogram, user selectable
- **Pulse Strength (Optional)**
8-segment bar graph indicates logarithmic strength of patient's pulse. Not proportional to pulse volume.

CAPNOGRAPH

- **Measurement**
Non-dispersive IR absorption
- **Calibration**
Manual 2 point
- **Measurement Range**
0 to 10% CO₂ STPD (standard temperature and pressure dry)
- **Display Range**
0 to 100 mmHg; 0-13.3 kPa; 0-10% CO₂
- **Accuracy**
±2 mmHg or 4% (whichever is greater) (for <50 breaths/min)
- **System Response Time**
2.46s to 90% with 8 ft. patient line
- **N₂O Compensation**
Selectable 40% (default = OFF)
- **Averaging**
4-breath average
- **Flow Rate**
120 ± 20 ml/min

RESPIRATION RATE

- **Range**
0 to 150 breaths/min
- **Accuracy**
±1 bpm
- **Averaging**
4-breath average

SpO₂ (OPTIONAL)

- **Range**
0 to 100% SpO₂ (functional)
- **Accuracy**
±2% at 70 to 100% SpO₂
±3% at 50 to 69% SpO₂
- **Averaging**
8 beats

PULSE RATE (OPTIONAL)

- **Range**
30 to 254 bpm
- **Accuracy**
±2 bpm or 2% (whichever is greater)
- **Averaging**
8 beats

ALARMS

- **Physiological**
User-selectable high and low limits, adjustable volume, screen indicators
- **System**
Adjustable volume
- **Silence**
Temporary 2 min or indefinite (permanent) LED indicator

SERIAL OUTPUT

Infrared Port compatible with MCP8850B Printer or similar device.

- **Data Format**
Non-standard format - 500 baud
- **Options**
Text only, no graphics. User-selectable patient data log or trend tables.

POWER REQUIREMENTS

- **AC Power**
120 V~, 60 Hz; 230 V~, 50/60 Hz
100 V~, 60 Hz (optional)
- **Battery Charger**
MoliEnergy

BATTERY LIFE

Li+ (lithium-Ion), 7.4 VDC.
User replaceable, internal rechargeable. Fully charged continuous life of approximately 6 hours. Maximum full-capacity charging is 2.5 hours.

DIMENSIONS

- **Base Unit**
5.0" H x 4.38" W x 2.9" D
(12.2 cm H x 11.1 cm W x 7.37 cm D)

WEIGHT

22 oz. (0.63 kg)

ENVIRONMENTAL

- **Operating Temperature**
32°F to 122°F (0°C to 50°C)
- **Storage Temperature**
-40°F to 167°F (-40° to 75°C)
- **Relative Humidity** (non-condensing)
15-95% (storage)
10-95% (operating)

STANDARDS

EN60601-1 (IEC 60601-1), UL544, EN60601-1-2 (IEC 60601-1-2)

REGULATORY COMPLIANCE

ISO 9001/ISO 13485

CATALOG NUMBER

- **8400**
Capnocheck® II Hand-Held Capnograph/Oximeter
- **8401**
Capnocheck® II Hand-Held Capnograph

Information subject to change without notice. Please read the Instructions for Use supplied with the product for detailed instructions, warnings and cautions.

Q-Line for U.S.A. customers: call 1-800-348-1842 for clinical application information.

For more information, please call Smiths Medical PM, Inc., at 262-542-3100 or 800-558-2345, or your Smiths Medical distributor.

Smiths Medical PM, Inc.

Patient Monitoring and Ventilation

N7W22025 Johnson Drive, Waukesha, Wisconsin 53186 USA

Phone: 262-542-3100 Fax: 262-542-0718 Toll-Free: 800-558-2345

www.smiths-medical.com/bci

Smiths Medical, part of the global technology business Smiths Group

BCI, Capnocheck and the Smiths Medical design mark are trademarks of the Smiths Medical family of companies. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trade names, trademarks or service marks of their respective owners. ©2008 Smiths Medical family of companies. All rights reserved. BCI 8400 Rev.06 05/08

smiths medical