

**NEW!  
KEYPAD FEATURE!**

# LifeDop™ 300 ABI


by SummitDoppler

## Economical Hand-Held Doppler System with FASTER ABI Calculation & Waveform Printout



Perform the ABI exam for P.A.D. diagnosis  
FASTER & EASIER with the LifeDop 300 ABI!

- Keypad for easily entering pressure data & calculating the indices
- Stores pressure data & two ankle waveforms prior to printing
- Proven technology: bi-directional probe and high speed printer
- Improved flow direction indicator on display

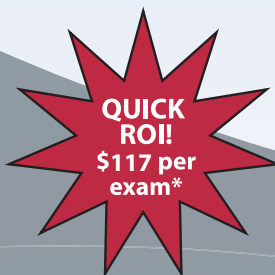

**To order the Life Dop 300 ABI use:**

**Model # L300AC: ABI System with 4 Cuffs & Aneroid**

\* The ABI exam for the diagnosis of peripheral arterial disease (P.A.D.) is reimbursable under CPT code 93922 at an average of \$117/exam based on published Medicare Fee Schedules, Part B Contractors, 2010. Arterial exams must be medically necessary in order to seek reimbursement. Other requirements, such as vascular accreditation may apply. Healthcare providers are advised to contact their third-party payers for specific coverage, coding and payment information.

**Summit Doppler Systems, Inc.**  
4680 Table Mountain Dr. # 150  
Golden, Colorado 80403  
#800-554-5090 or #303-423-7572  
Fax #303-940-7165  
[www.summitdoppler.com](http://www.summitdoppler.com)

# LifeDop™ 300 ABI

## Economical Hand-Held Doppler System with FASTER ABI Calculation & Waveform Printout



### ABI 4 Cuff Package includes:

- Doppler main unit
- 8 MHz bi-directional probe
- Portable printer with cable
- Aneroid
- (2) 10 cm cuffs
- (2) 12 cm cuffs
- Adhesive-backed label paper
- Report forms
- ABI training video
- ABI chart
- Ultrasound gel (60 g)
- Carrying case

SummitDoppler

**QUICK  
ROI!  
\$117 per  
exam\***

**To order the LifeDop 300 ABI use:**

**Model # L300AC: ABI System with 4 Cuffs & Aneroid**

### Main Unit and Probe

- Bi-directional 8 MHz high sensitivity probe (optional 5 MHz probe)
- Separate waveforms indicate flow toward and away from the probe
- Auto-scaling and no calibration required
- 200-4,000 Hz bandwidth covers peripheral applications
- User replaceable probe cable
- Low battery indicator
- Long battery life (3 AA alkaline batteries)
- Dimensions: 147 x 70 x 35 mm (h w l)
- Weight: 320 grams

### ABI Calculation and Waveform Storage

- Calculates ABI for both sides using the higher of the two arm pressures
- Accepts pressures from 30 to 299 mmHg
- Stores Doppler waveform from each ankle for printing

### Compact Portable Printer

- Label paper for fast report generation
- High speed
- Rechargeable (Li-Ion battery)
- Size: 136 x 84 x 43 mm (h w l)
- Weight: 280 grams



**Perform the ABI exam for P.A.D. diagnosis  
FASTER & EASIER with the LifeDop 300ABI!**

**Summit Doppler Systems, Inc.**

**4680 Table Mountain Dr. # 150**

**Golden, Colorado 80403**

**#800-554-5090 or #303-423-7572**

**Fax #303-940-7165**

**www.summitdoppler.com**

*\* The ABI exam for the diagnosis of peripheral arterial disease (P.A.D.) is reimbursable under CPT code 93922 at an average of \$117/exam based on published Medicare Fee Schedules, Part B Contractors, 2010. Arterial exams must be medically necessary in order to seek reimbursement. Other requirements, such as vascular accreditation may apply. Healthcare providers are advised to contact their third-party payers for specific coverage, coding and payment information.*

MKT0183B