


**New** Single Patient  
Surface Electrode

**New** Single Patient  
Surface Electrode

## Ambu® Neuroline 715

**Solid Gel**

**715**

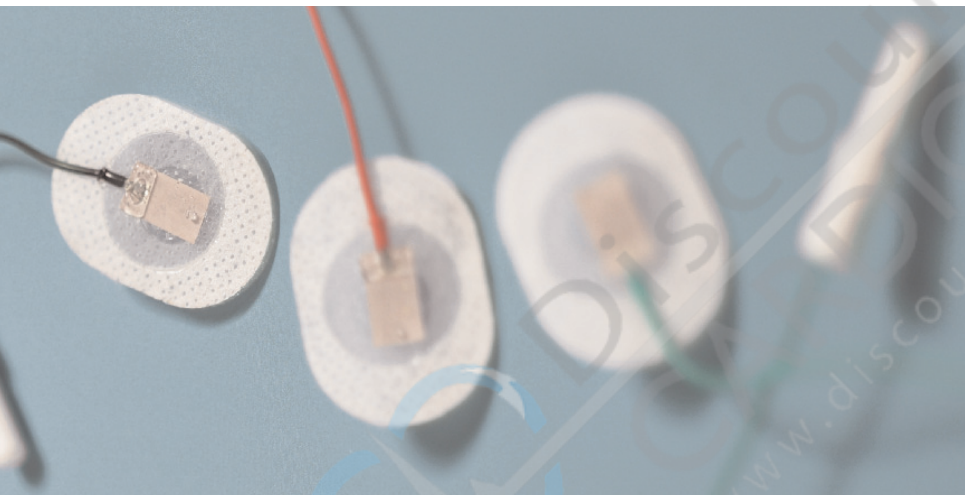

- Electrode size - 474 mm<sup>2</sup>
- Sensor area - 54 mm<sup>2</sup>

### The Easy and Effective Way to Record Muscle and Nerve Activity

The Ambu Neuroline 715 electrode features a high quality silver/silver chloride sensor for optimum signal quality. The transparent sensor area enables you to locate the examination target area very precisely. Moreover, the electrode comprises a skin friendly

solid-gel that allows you to reposition the electrode several times without losing its adhesive power. The Ambu Neuroline 715 electrode is conveniently available with colour-coded lead wires in a green/red/black combination and offered with different choices of lead wire lengths in lengths of 50, 80 or 120 cm.

### Unique features/benefits

- High quality Ag/AgCl sensor gives you an optimum signal quality
- Transparent sensor area allows precise sensor placement over the target area and center of electrode is marked for measuring distance
- Skin friendly solid gel enables you to reposition the electrode several times
- Color-coded lead wires give you easy placement identification
- Wound lead wire gives you an easy handling where the lead wires do not tangle

