



Key Features

- Proprietary design reduces pressure points on the forehead and across the bridge of the nose
- Built-in chin flap eliminates need for chin strap
- Less claustrophobic design improves patient comfort and eliminates nasal bridge pressure sores
- 3 sizes of pillows included to fit all patients' nostrils
- Easy to adjust and quick to fit headgear
- Reduces waste by identifying proper size before placement
- 100% latex free

Hybrid AV and NE

Hybrid disposable dual airway mask

Delivers optimal performance. Maximizes fit and comfort

The Hybrid AV & NE interfaces are disposable dual airway leak-free masks, designed for providing pressure support therapy or non-invasive ventilation. The unique design of the Hybrid promotes patient comfort.

The Hybrid AV features an anti-asphyxiation valve that allows for access to room air if pressure decreases while the Hybrid NE eliminates it, making it suitable for use with ventilators equipped with an active safety valve.

Each Hybrid package incorporates 3 sizes of pillows and a built-in mouth sizing gauge, eliminating the hassle of opening multiple packages to get the right fit.

SIZES

Small
Medium
Large

MATERIAL

Cushion	Silicone
Pillows	Silicone
Valve	Silicone
Caps	Silicone
Shell	Polycarbonate
Elbow	Polycarbonate
Rings	Polycarbonate
Buckles	Polypropylene
Headgear	Nylon, Neoprene, Velcro Hook, Loop




ORDERING INFORMATION

RespCare nr.	Model	Size	Order nr.
AV – for use with bilevel ventilatory support			
HYB551	Hybrid AV Disposable Mask	Small	417 003 000
HYB553	Hybrid AV Disposable Mask	Medium	417 003 001
HYB555	Hybrid AV Disposable Mask	Large	417 003 002



NE – for use with NIV mode of ICU ventilator

HYB571	Hybrid NE Disposable Mask	Small	417 004 000
HYB573	Hybrid NE Disposable Mask	Medium	417 004 001
HYB575	Hybrid NE Disposable Mask	Large	417 004 002

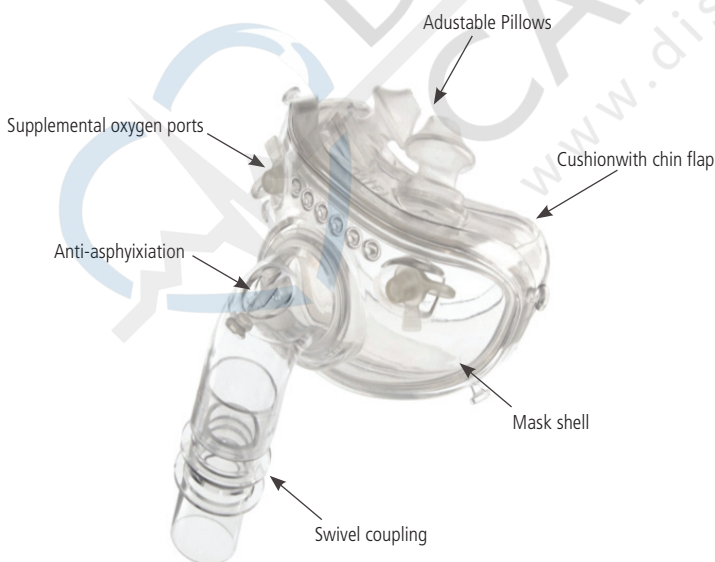

Manufactured for

RESCARE



US: Rx only

Distributed by

Ambu A/S

Baltorpbakken 13
DK-2750 Ballerup
Danmark

Tlf.: +45 72 25 20 00

Fax: +45 72 25 20 53

www.ambu.com