



OUR NEWEST MUSCLE STIMULATOR DELIVERS WORLD-CLASS PERFORMANCE.

Introducing the muscle stimulator that leads the field in every event. The Sys*Stim® 226 meets therapy challenges with ease—especially in clinics with wide-ranging needs.

Flexibility is key to the Sys*Stim 226's winning ways. Treatment waveforms include Interferential, Pre-modulation, Medium Frequency and Biphasic with adjustable duration. For added versatility, choose from Surge, Reciprocation, Surge with Reciprocation and Vector Rotation modes. What's more, the Sys*Stim 226 offers combination therapy with any waveform, giving you the ability to control independent therapies on different parts of the body, or even on different patients.

Whether at home or on the international circuit, the Sys*Stim 226 transcends language barriers with its user-friendly, WYSIWYG waveform graphics. It also conforms to worldwide safety requirements, so it can benefit your patients wherever you take it—in its small, lightweight tote.

At the finish line, the full-featured Sys*Stim 226 is a hands-down winner in value. So call (800) 854-9305. And see what a world-class performer can do for *your* clinic.



METTLER ELECTRONICS®

©2000, Mettler Electronics Corp. All rights reserved.

THE SYS*STIM 226 MUSCLE STIMULATOR




Numeric Keypad allows quick entry of time, frequency and other treatment protocols.



Select the **Treatment Waveform** by pressing the picture of the waveform button.



Tactile Membrane Entry Panel simplifies treatment control.

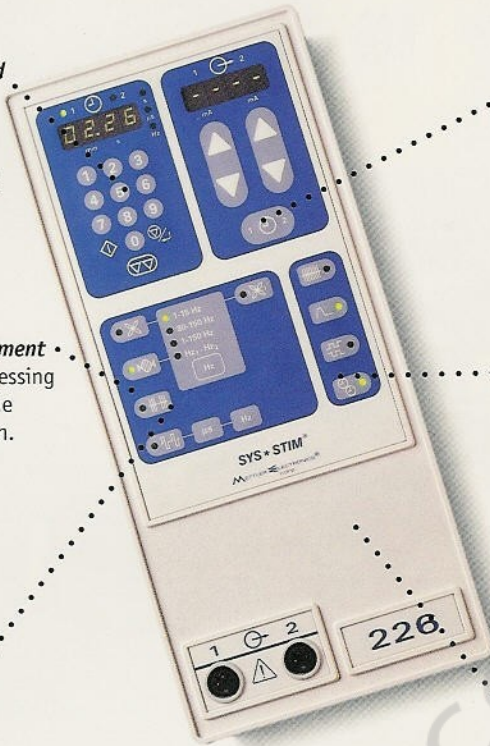

Program **On/Off** times and amplitude modulation with a single touch of a button.




Control two independent treatment protocols using the **Two Timer Mode**.




Surface Mount Technology and **Microprocessor control** improves reliability.




Combine any Sonicator® therapeutic ultrasound with the Sys*Stim 226 to perform combination therapy.

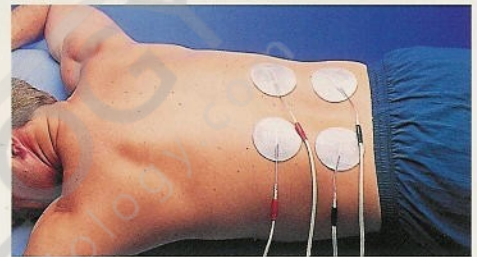

Treat low back pain using the interferential mode with amplitude modulation.

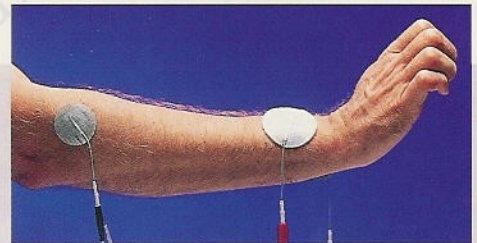

Re-educate two different muscle groups by combining the surge and reciprocation modes.

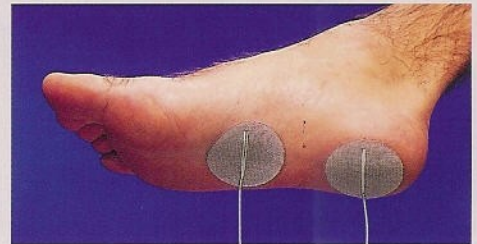

The pre-modulated waveform may even be used to stimulate the bottom of the foot.



The lightweight and portable Sys*Stim 226 can be easily taken to the patient.

Specifications

Input	115 VAC ± 10%, 50/60 Hz, 0.4 Amp Nom.
Weight	4.8 pounds
Dimensions	4.3 in (H) x 6 in (D) x 13.4 in (L)
ETL and C-ETL Listed	Model ME 226 (9801427)
Warranty	2 year limited
Waveforms	Interferential, pre-modulated, medium frequency (Russian), biphasic
Output	Interferential, 0-65 mA RMS, 1 Kohm load Pre-modulated, 0-50 mA RMS, 1 Kohm load Medium frequency, 0-50 mA RMS, 1 Kohm load Biphasic, 0-99 mA peak, 1 Kohm load
Amplitude Modulation	With interferential ±50% (Vector rotation) Surge, On/Off times programmable, Reciprocation, Equal On/Off times or use with Surge
Frequency	Interferential, 4000 Hz, 0-250 bps Pre-modulated, 1-250 bps Medium frequency, 2500 Hz, 50 pps Biphasic, 1-120 pps

Mettler Electronics Corp., 1333 S. Claudina St. Anaheim CA 92805
(800) 854-9305 Department J Fax (714) 635-7539
Tel. (714) 533-2221 <http://www.mettlerelec.com>

© Copyright 2000, Mettler Electronics Corp. All rights reserved.
Sonicator and Sys*Stim are registered trademarks of Mettler Electronics Corp.

METTLER ELECTRONICS®



IR7-86 Printed in U.S.A.
10M-10/00