



PageWriter 100, 200/300pi Cardiograph Series

Superior ECGs

by any measure

Let's make things better.



PHILIPS


Reliable, economical performance anywhere you need it




The Philips PageWriter 100 gives you the portability you need without compromising the features and functions essential for high-quality ECG recordings.



Everything about the PageWriter 100 makes ECG recording easy, including its patient cable which reduces lead tangling and quickens patient prep.

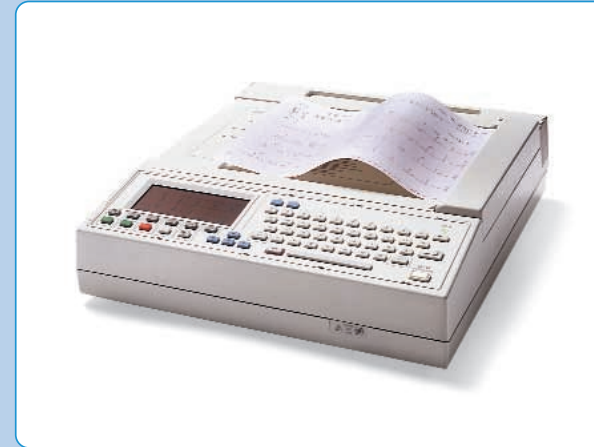

The alphanumeric keyboard on the Philips PageWriter 200/300pi lets you enter patient IDs in seconds while the display lets you check signal quality prior to printing, saving time and paper.

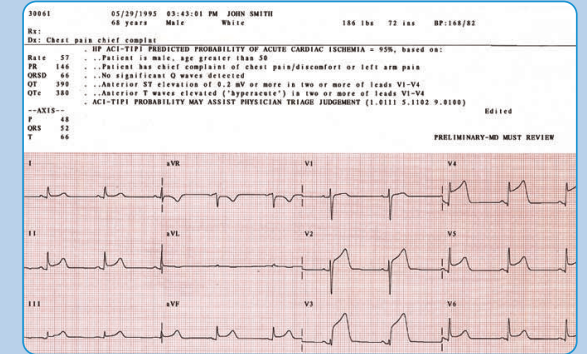

The PageWriter 300pi offers automatic interpretation of ECG data and advanced predictive instruments for fast patient assessment.

PHILIPS' PAGEWRITER 100 AND 200/300PI

For superior performance anywhere you need it—office, clinic or hospital—Philips offers high-quality cardiographs that range in function and price, so you can choose the model that best meets your needs and fits your budget.

The Philips PageWriter 100 is a cardiograph packed with features to make the most of your investment. Or move up in capability for an unbeatable combination of quality and value with the Philips PageWriter 200/300pi.

Either way you're assured of Philips' legendary reliability and low cost-of-ownership.

EASY ECONOMICS

The Philips PageWriter 100, our most economical cardiograph, is so easy to use there's virtually no learning curve. Activate an ECG at the press of a

single button and be assured of accurate results. Traces are printed on full-size paper from a continuous-feed, high-resolution, digital-array printer, providing maximum quality output—with no cutting or pasting required.

The PageWriter 100 records 3 or 6 channels at speeds ranging from 5 mm to 50 mm per second, so you can capture multiple levels of detail. Full diagnostic frequency meets American Heart Association (AHA) guidelines. Dedicated function keys identify frequently used operations in easy-to-understand terms. And drop-in paper loading eliminates the hassle of threading.

Powered by battery or AC, the PageWriter 100 is an ideal solution for those who need mobility in recording ECGs.

ADVANCED FEATURES, UNBEATABLE VALUE

The Philips PageWriter 200/300pi offers a wealth of advanced features designed to save you time and simplify the acquisition, interpretation, storage and transmission of ECGs.

As with the PageWriter 100, this standard cardiograph captures accurate 12-lead ECGs on full-size paper—with no cutting or pasting. Plus, it features an alphanumeric keyboard so you can enter patient IDs in seconds.

Move up to the PageWriter 300pi to get the same highly accurate adult and pediatric interpretive algorithm that's used in the Philips PageWriter Xli cardiograph. And the 300pi has a built-in ECG Analysis Program to produce reports that include ECG measurements and analysis statements.

The PageWriter 300pi also offers advanced predictive instrument programs—especially helpful to physicians making critical decisions in Emergency Departments.

ACI-TIPI (Acute Cardiac Ischemia-Time Insensitive Predictive Instrument) helps you predict the probability of acute cardiac ischemia, while TPI (Thrombolytic Predictive Instrument) helps you identify patients having acute myocardial infarction who might benefit from thrombolytic therapy.

Full-featured and lightweight, the PageWriter 200/300pi offers an unbeatable combination of quality and value for offices, clinics and hospitals.

Philips Medical Systems is part of Royal Philips Electronics

United States
Philips Medical Systems
Cardiac and Monitoring Systems
3000 Minuteman Road
Andover, MA 01810
(800) 934-7372

Canada
Philips Medical Systems
Canada
281 Hillmount Road
Markham, ON
L6C 2S3
(800) 291-6743

Europe, Middle East and Africa
Philips Medizinsysteme
Böblingen GmbH
Cardiac and Monitoring Systems
Hewlett-Packard Str. 2
71034 Böblingen
Germany
Fax: (+49) 7031 463 1552

Latin America Headquarters
Philips Medical Systems
1550 Sawgrass Corporate Parkway #300
Sunrise, FL 33323
Tel: (954) 835-2600
Fax: (954) 835-2626

Asia Pacific Headquarters
Philips Medical Systems
30/F Hopewell Centre
17 Kennedy Road
Wanchai, Hong Kong
Tel: (852) 2821 5888
Fax: (852) 2527 6727

www.medical.philips.com

medical@philips.com

fax: +31 40 27 64 887

Philips Medical Systems
Global Information Center
I.B.R.S. / C.C.R.I.
Numéro 11088
5600 VC Eindhoven
Pays-Bas /
The Netherlands
(no stamp required)

© 2002 Philips Electronics
North America
Corporation. All rights
are reserved. Reproduction
in whole or in part is
prohibited without the prior
written consent of the
copyright holder.

Printed in the United States.
June 1, 2002
5990-0522EN



HIGH QUALITY WITHOUT THE HIGH PRICE

Whether you choose the PageWriter 100 or 200/300pi, you will get the high-quality operation you expect from a Philips solution:

Outstanding ECG Quality

- Simultaneous acquisition of up to 12 leads
- Recording speeds of 5, 10, 25 or 50 mm/sec for capturing multiple levels of detail

Greater Flexibility

- Full-page English (A size) or metric (A4 size) reports for convenient reading and filing
- Versatile report formats and speed options quickly provide auto reports or rhythm reports
- Convenient battery operation for greater mobility
- Copy button lets you generate as many ECG copies as you want for filing and distribution
- Lightweight enough to move easily around office or clinic

Ease of use

- Color-coded keys make operation intuitive
- Easy loading allows fast, hassle-free paper replacement
- Provides money-saving alternative to time-consuming cutting, pasting and mounting of single-channel strips
- Dedicated function keys for frequently used operations include Chart, Speed, Size, V Leads, Filter, Page Advance, Auto, Manual, Copy, Stop — and on the 200/300pi there's a dedicated key for ID.

THE PHILIPS PAGEWRITER 100

Your investment in the PageWriter 100 also provides:

- Simple, one-button operation
- Continuous rhythm-strip recording of 3 or 6 channels
- Rechargeable lead acid battery to record approximately 40 auto ECGs and 40 copies or 40 minutes of rhythm
- An easy-to-clean system

THE PHILIPS PAGEWRITER 200/300PI

Choose the PageWriter 200/300pi to get:

- Real-time, continuous recording of 3, 6 or 12 channels
- Extensive ECG quality control features, including AC noise filtering and baseline wander correction
- User-configurable filters
- A keyboard to quickly and easily add patient ID information
- A display to preview the signal prior to printing, saving time and paper

ADVANCED INTERPRETATION AND ANALYSIS (PAGEWRITER 300PI)

When you choose the PageWriter 300pi you get automatic interpretation of ECG data and advanced predictive instruments for fast patient assessment:

- ACI-TIPI helps predict the probability of acute cardiac ischemia
- TPI helps identify patients having acute myocardial infarction (AMI) who might benefit from thrombolytic therapy
- Automatic interpretation of ECG data using clinically proven algorithm
- Adult and pediatric analysis programs are standard

U.S. WARRANTY

A three-year cardiograph warranty with overnight replacement is standard on the Philips PageWriter 100 and 200/300pi.



PHILIPS